

100% MADE IN ITALY



EXTROSTYLE CATALOG 2025 / 26

extrostyle
medical shoes

WELCOME TO

EXTROSTYLE

ExtroStyle is a company in constant evolution, deeply rooted in the values of Italian craftsmanship and projected toward innovation and sustainability. For over 35 years, we have worked with passion to improve the quality of life for those who choose to move each day with energy, safety, and awareness.

Well-being is the guiding principle behind everything we do: from design to production, we place the individual and their need for comfort, protection, and performance always at the center. Our footwear is designed to safeguard sensitive feet during physical activity, minimizing joint stress and enhancing stability in every situation.

Every model is thoughtfully designed to deliver deep, lasting comfort that accompanies the wearer throughout the day, contributing to their

overall well-being and to the joy of movement. Because for us, **walking is living.**

In our factories located in Umbria and the Marche region, a fully Italian production chain combines meticulous attention to detail, artisanal expertise, and a constant commitment to research and development. Each step of the process is carefully monitored to ensure maximum quality, meeting the highest standards of safety and reliability.

CERTIFICATIONS & RELIABILITY

All ExtroStyle footwear has the CE mark and is certified as Class I Medical Devices (in compliance with MDR 745/17), registered at the Italian Ministry of Health. Our products also meet the requirements of UNI EN ISO 9001 and UNI EN ISO 13485 standards.



100% MADE IN ITALY





OUR DISTRIBUTORS

Our vision extends beyond national borders: today, ExtroStyle is present in over 30 countries worldwide, bringing the excellence of Italian craftsmanship to the service of well-being everywhere.



INDEX

THE PRESCRIPTION OF ORTHOSES BY RISK PHASE	4
INSOLE SYSTEM	5
ACTIVE LINE	7
ENERGY LINE	23
ENERGY 10	34
HORIZON LINE	39
COMFORT LINE	45
REHABILITATION LINE	62
BAMBOO	73
XXL	77
BRACES FOR DIABETIC FOOT AND POST-OPERATIVE LINE	79
EXTREME LINE	91

100% MADE IN ITALY



THE PRESCRIPTION OF ORTHOSES BY RISK PHASE

LEVEL OF DIABETIC RISK	PROBLEM	GUIDELINE FOR THE ORTHOTIC SOLUTION FOOTWEAR AND ORTHOPEDIC INSOLES		IDEAL PRODUCT
0 LOW RISK	Diabetic patient with normal conformed foot. No neuropathy.	Calibrated Stocking		Active Line Energy Line Comfort Line Rehabilitation Line
		FOOTWEAR Comfort	ORTHOPEDIC INSOLE No one or protective	
1 MEDIUM RISK	Neuropathy with light deformities.	Perfect fitting		Active Line Energy Line Rehabilitation Line
		FOOTWEAR Theraupedic to avoid stress caused by an inappropriate upper OUTSOLE Flexible bio-mechanic	ORTHOPEDIC INSOLE Protective and customized	
2 HIGH RISK	Neuropathy or arterial disease. Relevant deformities but without ulceration and/or previous amputation.	Accomodate deformities and rebalance pressure peaks		Active Line Energy Line Rehabilitation Line
		FOOTWEAR Theraupedic OUTSOLE Flexible bio-mechanic or semi-rigid	ORTHOPEDIC INSOLE Customized	
3 VERY HIGH RISK	Neuropathy or arterial disease. Deformities, ulceration and/or previous amputation.	Avoid reulceration. Accomodate deformities and rebalance pressure peaks		Active Line Energy Line Diabetic foot and post operative device
		FOOTWEAR Theraupedic or customized INSOLE Bio-mechanic rigid or semi-rigid customized orthopedic braces	ORTHOPEDIC INSOLE Customized	

EXTROSTYLE INSOLE SYSTEM

The insoles used for the ACTIVE and ENERGY lines are composed of a high-density open-cell foam that ensures **absorption and deabsorption** capacity, keeping the foot ventilated, dry and odourless, as well as ensuring a pleasant sensation when walking. The **ergonomic design** of the footbed offers high comfort, minimising the effort and tensions in the plantar support, both during walking and in a static position, dampening the impact of pressure during walking.



CHARACTERISTICS OF THE ENERGY LINE INSOLES

→3 mm footbed

To be used immediately after an operation on the forefoot in the presence of very large dressings and bandages.

→6 mm footbed

To be used during rehabilitation and in the presence of bandages and swellings of the forefoot.

→9 mm footbed

To be used as a standard footbed already inside the footwear.



FIBRE INSOLE SYSTEM

Fibre insoles and TECNICA footwear are great allies against diabetic foot issues, as they allow the authorised healthcare professional to use this footwear at all stages of risk, respecting all international protocols for the treatment and care of diabetic foot issues.

CHARACTERISTICS OF FIBRE INSOLES

S-Rigid insole

Extra-rigid fibreglass insole. It makes the outsole completely rigid, limiting the joint stress of the foot during walking and reducing significantly the pressure peaks on the forefoot, redistributing them correctly.

S-Semirigid insole

Semi-rigid fibreglass insole. It makes the outsole semi-rigid, supporting the foot without forcing it. The pressure peaks on the forefoot are reduced and the load on the joints is limited.




S-CLOUD INSOLE SYSTEM

The S-Cloud insole is supplied as an accessory for the ENERGY, REHABILITATION and Diabetic foot and post operative device lines, and for the Tecnica Stelt and Tecnica Biomed products. The innovative S-Cloud insole allows to remove little squares in the lower part of the entire insole surface thanks to the pre-incision of the 6 mm Evaform. The orthopaedic technician will identify the pieces to be removed to allow the

immediate discharge of pressures in the ulcerated area or in the area at risk of ulceration, allowing to create a precise zone of unloading of the necessary size in a precise point. The 3 mm self-modeling memory foam inside the insole, reduces the forces of the edges of the unloading area on the plantar area, distributing thus the body weight over a wider surface avoiding in this way pressure peaks.

CHARACTERISTICS OF S-CLOUD INSOLES

1. Breathable, antibacterial, antimicrobial, anti-odour and hypoallergenic On-Steam cover.
2. 3 mm self-moulding memory foam.
3. 10 mm x 10 mm modular inserts in Evaform, 6 mm thick.
4. Fibre adhesive base to stabilise the Evaform modules.

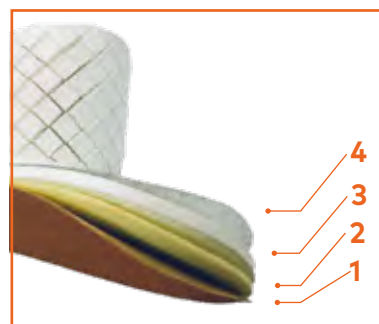

COMBINATION OF USE

The S-Cloud plantar system is divided into five lines of use:

- **S-Cloud Energy**
compatible with the ENERGY line (pair)
- **S-Cloud Fisio**
compatible with the REHABILITATION line (pair)
- **S-Cloud Diab**
compatible with the Diabetic foot braces line (piece)
- **S-Cloud Stelt** (piece)
- **S-Cloud Biomed** (piece)



S-AIR INSOLE SYSTEM

The S-Air insoles are included throughout the REHABILITATION LINE. They provide a pleasant

sensation when walking, avoiding also the consequences of the continuous impact on the joints.

CHARACTERISTICS OF S-AIR INSOLES

- Full ventilation thanks to the 100% breathable foam, which allows the air circulation.
- Keeps feet dry eliminating humidity through its absorption system.
- Eliminates odours, neutralising germs thanks to an antibacterial treatment.
- Prevents fatigue thanks to its shock absorption capacity.
- Longer-lasting foam thanks to its storage capacity.





ACTIVE LINE

ACTIVE LINE

INNOVATION AND PREVENTION FOR SENSITIVE FEET

The **ExtroStyle Active Line** was created to provide an advanced, practical solution for those who need specialized footwear for the prevention of diabetic foot and the management of common conditions such as hallux valgus, hammer toes, metatarsalgia, and rheumatic foot. It combines scientific research, functional design, and the meticulous craftsmanship that defines Italian quality.

The **Active Line**, which is the result of an advanced Research & Development process, is classified as a Class I Medical Device according to the European Directive MDR 745/17. It is developed in compliance with the latest international protocols in the field of diabetology.

These shoes represent the ideal combination of prevention and well-being, designed to protect sensitive feet, offer comfort during daily activities, and enhance the quality of life of the wearer.

Each model features a wide fit (12) and is designed without internal seams to prevent friction and microtrauma, helping to avoid the appearance of skin lesions. The materials used, which are breathable, self-molding, and strictly biocompatible, ensure an optimal fit and continuous protection, even during prolonged use.

The Active Line is designed to accommodate custom orthotics up to 6 mm, making the footwear adaptable to the clinical needs of each user. A distinctive technical feature is the presence of the Roller-profile outsole, designed with both anterior and posterior curvature and a flat central section extended to the metatarsals to improve balance. This design promotes the even distribution of the load, reducing pressure peaks on the forefoot and on the heel.

This structure helps to reduce the impact during walking, ensuring smooth, stable, and safe movement.



ACTIVE LINE WALKING MEANS LIVING.

The active **heel support** is padded, solid, envelops and protects the heel avoiding rubbing and injury.

The **PADDED LINING**, which is comfortable, breathable and **with no seams**, welcomes the foot in a protective and cushioning casing.

Designed to accommodate custom made insoles up to 6 mm thickness

PADDED LINING

Active outsole with roller profile to improve balance and promote even weight distribution

The central part of the outsole, being wide, guarantees **real stability** with consequent relaxation for the entire musculature.

The Active Line is available in THREE RIGIDITY CONFIGURATIONS – FLEXIBLE, SEMI-RIGID, AND RIGID – to specifically address different levels of preventive and functional needs. Whether for diabetic patients without lesions, individuals with foot deformities, or simply those seeking a reliable shoe for an active lifestyle, this line offers tailored solutions.

Thanks to its versatile use, careful material selection, and prevention-oriented design, the Active Line has become a benchmark for those seeking a technical and comfortable shoe, capable of protecting the foot and promoting everyday well-being.

BECAUSE PREVENTION IS LIVING, AND WALKING IS THE FIRST STEP TOWARD WELL-BEING!

PLANTARE S-RIGID ACTIVE

Certification no. 2702647

Sizes from 35 to 46

Ref. no. EX051710



PLANTARE S-SEMIRIGID ACTIVE

Certification no. 2702648

Sizes from 35 to 46

Ref. no. EX051710



PLANTARE ACTIVE 6 mm

Certification no. 2702649

Sizes from 35 to 46

Ref. no. EX050800



ACTIVE LINE FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDLED
LINING

ACTIVE MELISSA

Certification no. 2702629


Color Beige

Sizes from 35 to 43

Ref. no. EX057800

Diabetic foot risk level

0 1 2 3 see page 4

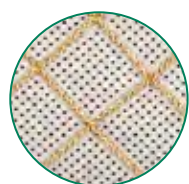

TOTAL OPENING
WITH DOUBLE ZIP

ACTIVE KATIA

Certification no. 2702630


Color Beige

Sizes from 35 to 43

Ref. no. EX057300

Diabetic foot risk level

0 1 2 3 see page 4

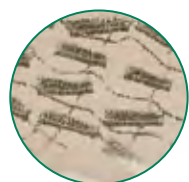

PADDLED AND
SOLID ACTIVE
HEEL SUPPORT

ACTIVE GORIZIA

Certification no. 2258592


Color White / Copper

Sizes from 35 to 43

Ref. no. EX057300

Diabetic foot risk level

0 1 2 3 see page 4

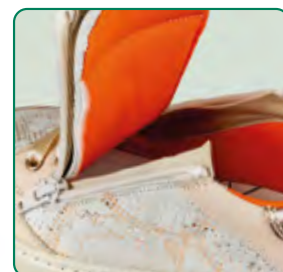

TOTAL OPENING
WITH DOUBLE ZIP

WOMAN

100% MADE IN ITALY

ACTIVE LINE

FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE NAXOS

Certification no. 2702631

Color Natural

Sizes from 35 to 43

Ref. no. EX057400

Diabetic foot risk level

0 1 2 3 see page 4



LATERAL OPENING WITH ZIP



UPPER IN FABRIC AND GENUINE
LEATHER



PADDED AND SOLID ACTIVE
HEEL SUPPORT

WOMAN

100% MADE IN ITALY

ACTIVE LINE

FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDLED
LINING

ACTIVE SOFIA

Certification no. 2755211

Color Black

Sizes from 35 to 43

Ref. no. EX057300

Diabetic foot risk level

0 1 2 3 see page 4



TOTAL OPENING
WITH DOUBLE ZIP



UPPER IN GENUINE LEATHER



PADDLED AND SOLID ACTIVE
HEEL SUPPORT



ACTIVE VIESTE PANNA DONNA

Certification no. 2420300

Color Cream

Sizes from 35 to 43

Ref. no. EX057450

Diabetic foot risk level

0 1 2 3 see page 4



TOTAL OPENING
WITH DOUBLE ZIP



ACTIVE LINE FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE LINE

FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE SCALEA NERO DONNA

Certification no. 2420301

Color Black

Sizes from 35 to 43

Ref. no. EX057700

Diabetic foot risk level

0 1 2 3 see page 4



100% MADE IN ITALY AND ATTENTION
TO EVERY DETAIL



TOTAL OPENING
WITH DOUBLE ZIP



EXTRA FLEXIBLE ANCHE SHOCK
ABSORBING ACTIVE OUTSOLE

WOMAN

100% MADE IN ITALY

ACTIVE LINE FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE MARATEA NERO DONNA

Certification no. 2420302

LITTLE CHAIN



Color Black

Sizes from 35 to 43

Ref. no. EX057950

Diabetic foot risk level

0 1 2 3 see page 4

ACTIVE ENNA NERO/ORO DONNA

Certification no. 2702646

GLITTER



Color Black/Gold
Glitter

Sizes from 35 to 43

Ref. no. EX057150

Diabetic foot risk level

0 1 2 3 see page 4

ACTIVE TROPEA NERO/ORO DONNA

Certification no. 2702644

GLITTER



Color Black/Gold
Glitter

Sizes from 35 to 43

Ref. no. EX057350

Diabetic foot risk level

0 1 2 3 see page 4

ACTIVE LINE

FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE NICOLAS

Certification no. 2702633

BREATHABLE UPPER



Color Natural

Sizes from 39 to 46

Ref. no. EX056900

Diabetic foot risk level

0 1 2 3 see page 4



LATERAL OPENING
WITH ZIP

ACTIVE DAVIDE

Certification no. 2761392


Color Green

Sizes from 39 to 46

Ref. no. EX057450

Diabetic foot risk level

0 1 2 3 see page 4

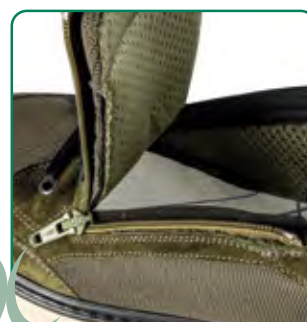

TOTAL OPENING
WITH DOUBLE ZIP



ACTIVE LINE
FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE DIEGO

Certification no. 2755213



Color Black

Sizes from 39 to 46

Ref. no. EX057450

Diabetic foot risk level

0 1 2 3 see page 4



TOTAL OPENING
WITH DOUBLE ZIP

ACTIVE LEONARDO

Certification no. 2755214



Color Taupe

Sizes from 39 to 46

Ref. no. EX057600

Diabetic foot risk level

0 1 2 3 see page 4



EXTRA FLEXIBLE
ANCHE SHOCK
ABSORBING ACTIVE
OUTSOLE

ACTIVE LINE
FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

**ACTIVE MILAZZO
GRIGIO UOMO**

Certification no. 2258596



Color Black/Grey

Sizes from 39 to 46

Ref. no. EX056900

Diabetic foot risk level

0 1 2 3 see page 4

**ACTIVE PALERMO
NERO UOMO**

Certification no. 2258590



Color Black/Bordeaux

Sizes from 39 to 46

Ref. no. EX056900

Diabetic foot risk level

0 1 2 3 see page 4

**ACTIVE ELIANO
BLU UOMO**

Certification no. 2702640



Color Blue/Grey

Sizes from 39 to 46

Ref. no. EX057350

Diabetic foot risk level

0 1 2 3 see page 4

ACTIVE LINE FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE STILO GRIGIO UOMO

Certification no. 2258604

Color Grey

Sizes from 39 to 46

Ref. no. EX057450

Diabetic foot risk level

0 1 2 3 see page 4



ACTIVE CATANIA NERO UOMO

Certification no. 2702636

Color Nero/T.Moro

Sizes da 39 a 46

Ref. no. EX057150

Diabetic foot risk level

0 1 2 3 see page 4



RIGID
OUTSOLE

ACTIVE MESSINA NERO UOMO

Certification no. 2702635

Color Nero/Dark brown

Sizes from 39 to 46

Ref. no. EX057900

Diabetic foot risk level

0 1 2 3 see page 4



ACTIVE TRAPANI 14 T.MORO UOMO

Certification no. 2702638

FIT
14

RIGID
OUTSOLE

Color Dark brown/
Blue/Brown

Sizes from 39 to 46

Ref. no. EX058150

Diabetic foot risk level

0 1 2 3 see page 4



MAN

100% MADE IN ITALY

ACTIVE LINE FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE ALCAMO NERO UOMO

Certification no. 2258597

Color Black/Grey

Sizes from 39 to 46


Ref. no. EX057350

Diabetic foot risk level

0 1 2 3 see page 4



BREATHABLE UPPER

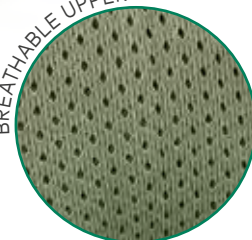

LACES FOR ADJUSTABLE OPENING



TOTAL OPENING
WITH DOUBLE ZIP



EXTRA FLEXIBLE ANCHE SHOCK
ABSORBING ACTIVE OUTSOLE

ACTIVE LINE
FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE ROCCO

Certification no. 2755212



Color Brown

Sizes from 39 to 46

Ref. no. EX057800

Diabetic foot risk level
0 1 2 3 see page 4



SIDE OPENING WITH ZIP AND
PADDED UPPER



GENUINE LEATHER



PADDED AND SOLID ACTIVE
HEEL SUPPORT

ACTIVE LINE FIT 12



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ACTIVE TARANTO NERO UOMO

Certification no. 2702641

Color Black/Bordeaux

Sizes from 39 to 46

Ref. no. EX057700

Diabetic foot risk level

0 1 2 3 see page 4



ACTIVE CONERO NERO UOMO

Certification no. 2258598

Color Black

Sizes from 39 to 46




Ref. no. EX057900

Diabetic foot risk level

0 1 2 3 see page 4



INTERNAL MEASURES

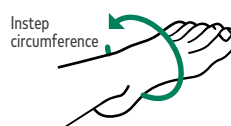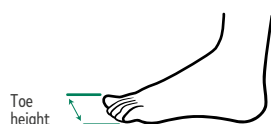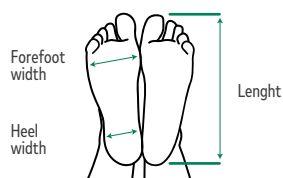

MEASURES OF THE MODELS: FORMIA, ALBA, SAVONA, GELA, VIESTE, CATANIA, CAPRI, PALERMO, MARATEA, SCALEA, ENNA, TROPEA, MILAZZO, ELIANO, STILO, ELIA, ALCAMO, TARANTO, CONERO.

FIT 12													
Sizes	UM	35	36	37	38	39	40	41	42	43	44	45	46
Lenght	mm	238	244	251	258	264	271	278	284	291	298	304	311
Forefoot width	mm	80	81,5	83	84	86	87,5	89	93	94	95,5	96,5	98
Heel width	mm	56	57	57,5	59	60	61	62	65	66	67	68	69
Forefoot circumference	mm	231	235,5	239,5	244	249,5	253	258	268	272	277	281	286
Toe height	mm	29	29,5	30	30,5	31	31,5	32	32,5	34	34,5	35	35,5
Instep circumference	mm	231	235,5	240	244,5	249	253,5	258	267,5	272	276,5	281	285,5

MEASURES OF THE MODELS: ACTIVE TRAPANI TESTA MORO UOMO

FIT 14													
Sizes	UM	35	36	37	38	39	40	41	42	43	44	45	46
Lenght	mm	238	244	251	258	264	271	278	284	291	298	304	311
Forefoot width	mm	83	84	85,5	87	88,5	90	92	95	97	99	100	102
Heel width	mm	57	58	59	60	61	62	63	66	67	69	69	70
Forefoot circumference	mm	241	245	251	255	260	264	268	280	284	291	294	298
Toe height	mm	31,5	32	32,5	33	33,5	34	34,5	37	37,5	37,5	38	39
Instep circumference	mm	241	245,5	250	254,5	259	263,5	268	280	284,5	289	293,5	298



ENERGY LINE

ENERGY LINE

ADVANCED THERAPEUTIC SUPPORT FOR DIABETIC FOOT MANAGEMENT



The Energy Line is ExtroStyle's most advanced solution for the therapeutic management of diabetic foot, designed to support patients through all risk levels, from grade 0 to grade 3. It also addresses other conditions such as hallux valgus, hammer toes, metatarsalgia, and rheumatic foot.

Classified as a Class I Medical Device under the European MDR 745/17 Regulation, it is a reliable reference for those in need of a technical, functional, and customizable shoe for the active care of pathological feet. Characterized by an extra-wide fit (fitting 14) and increased internal volume, the Energy Line is designed to accommodate feet with deformities, calluses, or bandages, assuring maximum comfort and reducing pressure points. The interior of the shoe is completely seamless to

eliminate any risk of friction or microtrauma: an essential feature in the prevention of ulcerations. Thanks to its design, which allows the insertion of customized orthotics up to 9 mm, each model can be adapted to the patient's specific needs, contributing to achieve an optimal plantar pressure distribution and significantly reducing the risk of complications.

A distinctive feature is the outsole with a roller profile, characterized by anterior and posterior curvature and a flat central section, which extends to the metatarsal area. This design promotes balance, pressure relief, and stability while walking. It enhances weight distribution, supports a smooth and safe gait, and reduces the risks associated with instability, pain, and ulcerations.

The roller-profile outsole is available in three different levels of rigidity — Flexible, Semi-Rigid, and Rigid — to perfectly meet the therapeutic needs of each patient:

- **ROLLER OUTSOLE ENERGY FLEXIBLE:** recommended for prevention (grades 0 and 1), in the absence of lesions, or for the treatment of rheumatic foot. It provides protection without restricting natural movement.
- **ROLLER OUTSOLE ENERGY SEMI-RIGID:** ideal for patients with risk level 2, experiencing pain, deformities, or previous ulcerations, who require greater stability and reduced pressure peaks.
- **ROLLER OUTSOLE ENERGY RIGID:** designed for high-risk patients (grade 3) with active or previous ulcerations, or in post-amputation phase. It ensures maximum pressure control and safe support while walking.

ENERGY IS MUCH MORE THAN A SHOE: IT IS AN INTEGRATED THERAPEUTIC TOOL DESIGNED TO PREVENT, PROTECT, AND TREAT DIABETIC FOOT AND OTHER PATHOLOGIES EFFECTIVELY, IN FULL COMPLIANCE WITH CLINICAL PROTOCOLS.

ENERGY LINE

CARE STARTS WITH THE RIGHT STEP

The reinforced, solid heel support surrounds and protects the heel, preventing rubbing and injuries.

The padded, comfortable, and breathable upper—designed without seams—envelops the foot in a protective and cushioning shell.

Designed to accommodate custom made insoles up to 9 mm thickness

Breathable upper

Roller-profile outsole with front and rear curvature

Improved weight distribution that promotes a smooth and safe gait

PLANTARE ENERGY

Certification no. 3 mm 2702686
6 mm 2702687
9 mm 2702688

Sizes from 35 to 50
Ref. no. EX030800

PLANTARE S-RIGID ENERGY

Certification no. 2702678

Sizes from 35 to 50
Ref. no. EX031710



PLANTARE S-CLOUD ENERGY

Certification no. 2702702

Sizes from 35 to 50
Ref. no. EX032220

PLANTARE S-SEMIRIGID ENERGY

Certification no. 2702701

Sizes from 35 to 50
Ref. no. EX031710



ENERGY LINE
FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING



ENERGY GIULIA

Certification no. 2760526

Color Beige

Sizes from 35 to 43

Ref. no. EX036800

Diabetic foot risk level

0 1 2 3 see page 4



ENERGY PETRA

Certification no. 2421897


Color Dove grey

Sizes from 35 to 43

Ref. no. EX036950

Diabetic foot risk level

0 1 2 3 see page 4

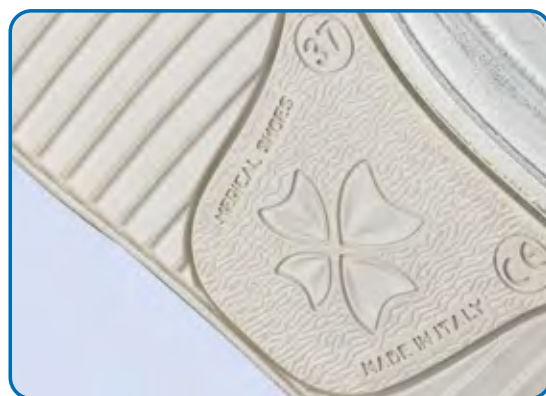

ENERGY LINE FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING



ENERGY ELOISA

Certification no. 2702662

Color Beige

Sizes from 35 to 43

Ref. no. EX037500

Diabetic foot risk level

0 1 2 3 see page 4



ENERGY LINDA

Certification no. 2702663

Color Black

Sizes from 35 to 43

Ref. no. EX037350

Diabetic foot risk level

0 1 2 3 see page 4



ENERGY LINE

FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ENERGY MARINA

Certification no. 2702692

Color Black

Sizes from 35 to 43

Ref. no. EX036900

Diabetic foot risk level

0 1 2 3 see page 4



ADJUSTABLE DOUBLE
VELCRO OPENING



ENERGY LINE

FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ENERGY ALICE

Certification no. 2702675



STRETCH
UPPER

Color Bronze

Sizes from 35 to 43

Ref. no. EX036550

Diabetic foot risk level

0 1 2 3 see page 4



COMFORTABLE,
HARMONIOUS AND ELEGANT

ENERGY ELENA

Certification no. 2702681



STRETCH
UPPER


Color Black

Sizes from 35 to 43

Ref. no. EX035800

Diabetic foot risk level

0 1 2 3 see page 4

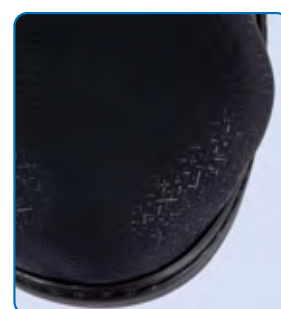

ENERGY LINE

FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ENERGY ALAN

Certification no. 2760527




Color Black

Sizes from 39 to 50

Ref. no. EX037800

Diabetic foot risk level

0 1 2 3 see page 4

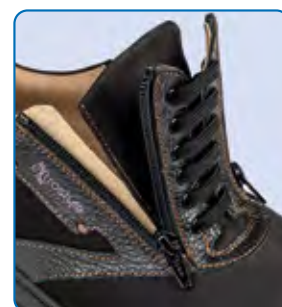

DOUBLE ZIP OPENING

ENERGY THOMAS

Certification no. 2702669



Color Black

Sizes from 39 to 50

Ref. no. EX036800

Diabetic foot risk level

0 1 2 3 see page 4



STRETCH
UPPER



VELCRO OPENING

**ENERGY IVAN**

Certification no. 2421894

Color Grey/Black

Sizes from 39 to 50

Ref. no. EX037550

Diabetic foot risk level

0 1 2 3 see page 4



REMOVABLE TECNICA INSOLE
TO COVER ALL RISK PHASES OF
THE SENSITIVE FOOT

**ENERGY LINE**
FIT 14REMOVABLE
INSOLEEXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICESUITABLE FOR THE
PREVENTION OF
DIABETIC FOOTPADDED
LINING

ENERGY LINE
FIT 14REMOVABLE
INSOLEEXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICESUITABLE FOR THE
PREVENTION OF
DIABETIC FOOTPADDED
LINING

ENERGY LORIS

Certification no. 2421682


Color Black

Sizes from 39 to 50

Ref. no. EX037100

Diabetic foot risk level

0 1 2 3 see page 4

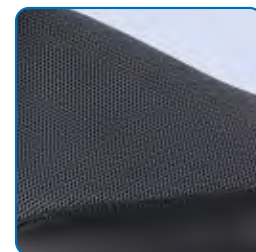

BREATHABLE UPPER

ENERGY LUIGI

Certification no. 2702694


Color Black

Sizes from 39 to 50

Ref. no. EX037100

Diabetic foot risk level

0 1 2 3 see page 4

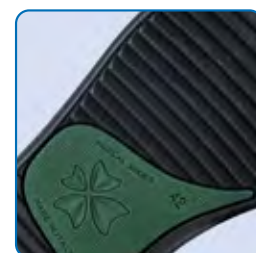STRETCH
UPPERBIOMECHANICAL
OUTSOLE

ENERGY VALENTINO NERO

Certification no. 2702693


Color Black

Sizes from 39 to 50

Ref. no. EX037900

Diabetic foot risk level

0 1 2 3 see page 4

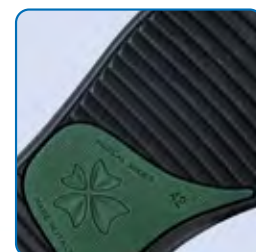BIOMECHANICAL
OUTSOLE

ENERGY LINE

FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ENERGY PAOLO MARRONE

Certification no. 2702700

Color brown
Sizes from 39 to 50
Ref. no. EX037100

Diabetic foot risk level
0 1 2 3 see page 4



ENERGY PAOLO NERO

Certification no. 2702684

Color Black
Sizes from 39 to 50
Ref. no. EX037100

Diabetic foot risk level
0 1 2 3 see page 4



ENERGY DANIELE

Certification no. 2702672



Color Black
Sizes from 35 to 50
Ref. no. EX038600

Diabetic foot risk level
0 1 2 3 see page 4



ENERGY ALEX NERO

Certification no. 2702691



Color Black
Sizes from 39 to 50
Ref. no. EX036800

Diabetic foot risk level
0 1 2 3 see page 4



ENERGY 10





PHASE 1 PREVENTION

Thanks to the predisposition for the Energy 10 **removable insole**, it is a great ally for the prevention or containment of the most common foot diseases such as: hallux valgus, hammertoes, metatarsalgia, heelitis and vascular diseases as well as for the prevention with a 0 - 1 risk in diabetic foot. The **flexible biomechanical energy outsole** allows a natural transition between the three phases of walking,

ensuring a reduction of the pressure peaks on the forefoot and limiting the wear on the joints, thus allowing a more spontaneous walk and reducing the feeling of fatigue. The **self-moulding and seamless upper** at the rubbing points envelops the anatomy of the foot and its possible deformities, ensuring comfort and preventing painful friction.



DIABETIC FOOT
VASCULOPATHIES



HALLUX VALGUS



HAMMER TOES
CLAW TOES
METATARSALGIA



RHEUMATIC FOOT
MORTON'S NEUROMA



BURSITIS
TARSAL TUNNEL



CHARCOT FOOT

PHASE 2 POST-OPERATIVE REHABILITATION

The energy 10 **full opening shoe with 2 zips**, equipped with a flexible outsole, can be defined as a multifunctional shoe thanks to a combination of insoles that are variable in height. **It can be worn from the first day after the operation until the end of convalescence.** This implies a more correct and balanced walking and distribution of body weight, as well as an undeniable benefit for the back, pelvis and legs. Also the psychological well-being is not to be underestimated, since the energy 10 footwear can be worn on both feet

and immediately after the operation. A **removable heel** (patented system) inclined by 8° can be applied to keep the foot in talism for all rehabilitation courses where the surgeon deems it necessary. The Energy 10 footwear was subjected by the Moticon Science (independent company) of Munich, to test the **PLANTAR DISCHARGE**, comparing it with five other footwear of other brands among the most sold on the European market and was the **BEST AND MOST EFFECTIVE IN ALL THE COMPARISON TESTS.**

TOTAL OPENING THANKS
TO THE 2 ZIPS

DIFFERENT HEIGHTS OF THE
INSOLES ACCORDING TO THE NEEDS

REMOVABLE HEEL PATENTED TO KEEP
THE FOOT IN TALISM



ENERGY 10, THE BEST EUROPEAN FOOTWEAR FOR PLANTAR DISCHARGE!

Study conducted by Moticon Science of Munich

PHASE 3 PROTECTION OF FOOT HEALTH

The Energy 10 footwear has been designed for men and women who want to live their natural "vital energy". The innovative technical characteristics open the way to a new idea of well-being, health and vitality. Indeed the line

is designed **to meet the needs and to prevent and manage the problems of the sensitive foot** but above all, to live well.

ENERGY LINE
FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING

ENERGY 10 BIO

Certification no. 2702664

FOOTWEAR
→PREVENTION
→POST-OPERATIVE
→ REHABILITATION
→FOR THE
PROTECTION OF
FOOT HEALTH



Color	Natural
Sizes	from 35 to 43
Ref. no.	EX036700

Diabetic foot risk level
0 1 2 3 see page 4


ENERGY 10



TOTAL OPENING WITH DOUBLE ZIP



SELF-MODELLING UPPER

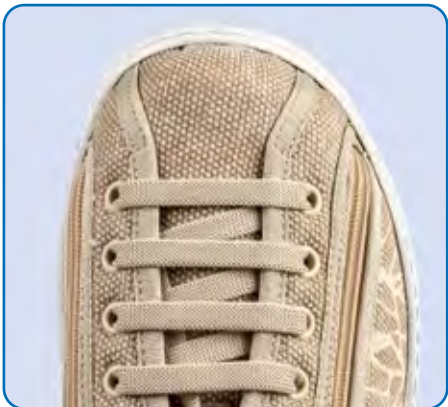

PUNTALE FILOFORMA

ENERGY LINE

FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



PADDED
LINING



TOTAL OPENING WITH
DOUBLE ZIP

ENERGY 10 BLU

Certification no. 2421679

Diabetic foot risk level

0 1 2 3 see page 4

Color Blue

Sizes from 35 to 41 (woman)
from 42 to 50 (man)

Ref. no. EX036150 (WOMAN)
EX036350 (MAN)



ENERGY 10 NERO

Certification no. 2702666

Diabetic foot risk level

0 1 2 3 see page 4

Color Black

Sizes from 35 to 41 (woman)
from 42 to 50 (man)

Ref. no. EX036150 (WOMAN)
EX036350 (MAN)



ENERGY 10 BIANCO

Certification no. 2702699

Diabetic foot risk level

0 1 2 3 see page 4

Color White / Black

Sizes from 35 to 41
(WOMAN)

Ref. no. EX036150



ENERGY PLUS

Certification no. 2702685

Color Black

Sizes from 35 to 46

Ref. no. EX031030



Sizes TABLE
35-36
37-38
39-40
41-42-43
44-45-46



ENERGY LINE

FIT 14

REMOVABLE
INSOLEEXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICESUITABLE FOR THE
PREVENTION OF
DIABETIC FOOTPADDED
LININGENERGY MORRIS
RIGID

Certification no. 2702695

RIGID
OUTSOLE

Color Blue

Sizes from 35 to 50

Ref. no. EX037700

Diabetic foot risk level

0 1 2 3 see page 4

ENERGY JACQUES
FLEX

Certification no. 2702698

Color Blue

Sizes from 35 to 50

Ref. no. EX035550

Diabetic foot risk level

0 1 2 3 see page 4






The **ENERGY MORRIS RIGID** and **ENERGY JACQUES FLEX**, thanks to the double velcro closure, have a complete opening in the anterior part of the shoe, making it more practical to wear. The extra wide fit and the adjustable internal volume make these shoes suitable for patients with large, bandaged, swollen or edematous feet. In addition, they are designed for the insertion of a customized insole.

→ **SEMI-RIGID OUTSOLE**: to be used in primary

prevention with risk degree 0 or 1, polyarthritis, Morton's neuroma, osteoarthritis, etc.

→ **RIGID OUTSOLE**: to be used in secondary prevention with risk degree 2 and 3, for patients with previous ulcerations at risk of recurrence, for patients with deformities, ulcerations or previous amputations, arthrodesis of the toes, arthrodesis of the tarsus, Charcot foot or biomechanical problems of movement.

INTERNAL MEASURES

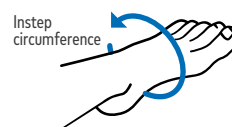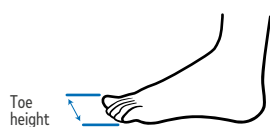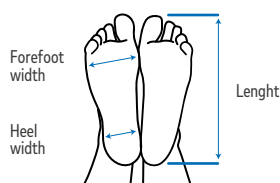

MEASURES OF THE MODELS: PERLA, AGATA, PETRA, VALENTINA, MARINA, ELENA, LUCIA, ALICE, IVAN, VALENTINO, LORIS, LUIGI, PAOLO, CARLO, DANIELE, ALEX, ENERGY 10, MORRIS, JACQUES.

Sizes	UM	35	36	37	38	39	40	41	42	43	44	45	46	47	48	49	50
Lenght	mm	233	240	247	253	260	267	273	280	287	293	300	307	313	320	327	333
Forefoot width	mm	84,5	86	87,5	89	90,5	92	93,5	96,5	98	99,5	101	102,5	104	105,5	107	108,5
Heel width	mm	56,4	57,6	58,8	60	61,5	62,4	63,6	66,8	68	69,2	70,4	71,6	72,8	74	75,2	76,4

INSOLE 3 mm																	
Forefoot circumference	mm	232,5	237	241,5	246,5	251	255,5	260	267	271,5	276	280,5	285	289,5	294	298,5	302,5
Toe height	mm	27	28	29	30	31	32	33	34	35	36	37	38	39	40	41	42

INSOLE 6 mm																	
Forefoot circumference	mm	227	231,5	236	240,5	245	249,5	254	261	265,5	270	274,5	279	283,5	288	292,5	297
Toe height	mm	24	25	26	27	28	29	30	31	32	33	34	35	36	37	38	39

INSOLE 9 mm																	
Forefoot circumference	mm	221	225,5	230	234,5	239	243,5	248	255	259,5	264	268,5	273	277,5	282	286,5	291
Toe height	mm	21	22	23	24	25	26	27	28	29	30	31	32	33	34	35	36



HORIZON LINE

HORIZON LINE

FUNCTIONAL ELEGANCE AND EVERYDAY WELL-BEING

The Horizon Line was created to combine comfort, prevention, and style, offering a refined and practical solution for those who wish to protect foot health without sacrificing elegance.

Ideal for the prevention of diabetic foot and for managing conditions such as hallux valgus, hammer toes, metatarsalgia, and rheumatic foot, this line is designed for everyday use.

The Horizon Line features a wide fit (fitting 11) and a widebase to ensure adequate space, even in the presence of deformities or orthotics. The interior is seamless, while the self-molding materials adapt

to the foot's shape, minimizing friction, localized pressure, and the risk of microtrauma. The Horizon outsole, with a 30 mm heel, is designed to ensure balanced weight distribution: body weight is evenly divided between the rearfoot and forefoot, promoting a smooth gait, correct posture, and lasting comfort. It is also slip-resistant to provide maximum safety on all types of surfaces.

Velcro or zip closures make the shoe easy to put on, even for individuals with limited mobility or difficulty using their hands.

THE ELEGANT AND DISCREET DESIGN MAKES THIS SHOE SUITABLE FOR DAILY USE IN ANY SETTING, FROM WORK TO LEISURE. **PREVENTION BEGINS WITH EVERYDAY WELL-BEING.**



**HORIZON CANAZEI**

Certification no. 2469201

Color Asphalt

Sizes from 35 to 43

Ref. no. EX047200

Diabetic foot risk level

0 1 2 3 see page 4



ELASTIC UPPER FOR DIFFERENT PATHOLOGIES

**HORIZON
LINE
FIT 11**REMOVABLE
INSOLEEXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICEPADDED
LINING

WOMAN

100% MADE IN ITALY

HORIZON LINE FIT 11



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

HORIZON ESTE

Certification no. 2469203

Color Brown

Sizes from 35 to 43

Ref. no. EX047550

Diabetic foot risk level

0 1 2 3 see page 4



TOTAL OPENING

HORIZON JESOLO

Certification no. 2469204

Color Black

Sizes from 35 to 43

Ref. no. EX047350

Diabetic foot risk level

0 1 2 3 see page 4



ZIP OPENING



HORIZON
LINE
FIT 11



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

HORIZON LEVICO

Certification no. 2469205

Color Sand

Sizes from 35 to 43

Ref. no. EX047300

Diabetic foot risk level

0 1 2 3 see page 4



STRETCH
UPPER



ADJUSTABLE ELASTIC
AND VELCRO OPENING

HORIZON CADORE

Certification no. 2469206

Color Brown

Sizes from 35 to 43

Ref. no. EX047500

Diabetic foot risk level

0 1 2 3 see page 4



ELASTIC UPPER NEAR
THE BIG TOE

HORIZON GARDA

Certification no. 2469207

Color Beige

Sizes from 35 to 43

Ref. no. EX047350

Diabetic foot risk level

0 1 2 3 see page 4



ELASTIC UPPER NEAR
THE BIG TOE

HORIZON LINE FIT 11



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING





MEDICAL
DEVICE



PADDED
LINING

INTERNAL MEASURES

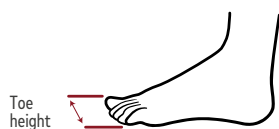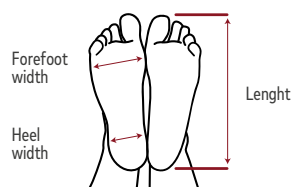

MEASURES OF THE MODEL: HORIZON CANAZEI - HORIZON CAORLE - HORIZON ESTE - HORIZON JESOLO - HORIZON CADORE - HORIZON GARDA.

Sizes	UM	35	36	37	38	39	40	41	42	43
Lenght	mm	234,6	241,3	248	254,6	261,3	267,9	274,6	281,3	287,9
Forefoot circumference	mm	223	227	232	236	241	245	250	254	259
Instep circumference	mm	225	229	234	238	243	247	252	256	261
Forefoot width	mm	81	82,5	84	85,5	87	88,5	90	91,5	93
Heel width	mm	52	53	54	55	56	57	58	59	60
Toe height	mm	18	19	20	21	22	23	24	25	26



COMFORT LINE

TECHNICAL AND STYLISH SANDALS AND SLIPPERS

Within the Comfort Line there is a collection of open footwear designed to provide maximum comfort and elegant style for both men and women. This line, which includes technical

and stylish sandals and slippers, is ideal for everyday use both indoors and outdoors, with special attention to the needs of those seeking functional footwear.

TECHNICAL FEATURES OF THE OPEN FOOTWEAR, COMFORT LINE

1. COMFORTABLE FIT

Line is designed to offer a wide and comfortable fit (fitting 10), able to accommodate different foot shapes, including sensitive feet or those with mild deformities.

2. ANATOMIC AND HIGH-PERFORMANCE INSOLE

The technical heart of these shoes is the anatomical insole, designed to provide optimal support:

- **Arch support:** Shaped to promote proper foot alignment and improve weight distribution.
- **Heel cup:** Stabilizes the foot and reduces impact during walking.
- **Metatarsal area:** Relieves pressure on the forefoot, ensuring comfort even during prolonged use.

This design makes the shoes ideal for preventing fatigue and improving the biomechanics of walking.

3. TECHNICAL AND SELF-MOLDING MATERIALS FOR THE UPPER

The upper is made of high-quality, self-molding, and breathable materials that adapt to the shape of the foot, ensuring personalized comfort. This feature is especially beneficial for those with deformities such as hallux valgus or for sensitive feet.

4. ELASTIC INSERTS ON THE UPPER

To provide additional comfort, elastic inserts are placed in the most sensitive areas (such as the big toe region) to reduce pressure and improve fitting, preventing irritation and discomfort.

5. DESIGNED FOR CUSTOM ORTHOTICS

All footwear of this line is designed to accommodate custom orthotics, providing a versatile solution that adapts to the specific needs of each customer.

6. SLIP-RESISTANT OUTSOLE

The outsole is made of high-traction materials, providing excellent safety on all surfaces. Additionally, it is lightweight but also durable, ideal for daily use, ensuring both stability and long-lasting performance.

THE OPEN FOOTWEAR OF THE COMFORT LINE REPRESENTS AN IDEAL COMBINATION OF TECHNICAL PERFORMANCE, ELEGANCE, AND COMFORT. WHETHER IT'S TECHNICAL SANDALS FOR OUTDOOR USE OR STYLISH SLIPPERS FOR HOME RELAXATION, THIS LINE OFFERS A VERSATILE AND HIGH-PERFORMING SOLUTION FOR THOSE SEEKING COMFORT, CUSTOMIZATION, AND STYLE—WITHOUT COMPROMISING FOOT WELL-BEING.

PLANTARE COMFORT

Certification no. 2473528

Sizes from 35 to 43

Ref. no. EX040880



PLANTARE CORK MEMORY

Certification no. 2473529

Sizes from 35 to 42

Ref. no. EX040950



PLANTARE CORK SLOW

Certification no. 2473530

Sizes from 39 to 50

Ref. no. EX040950



THE PLANTARE COMFORT CAN BE INSERTED IN: TECNICA ADELE PETROLIO, TECNICA ADELE ROSSO, TECNICA ADELE OCRA, TECNICA ADELE SUDAN, TECNICA ADELE FIORI, TECNICA ANNA FIORI, TECNICA SIMONA LILLA, TECNICA SIMONA VERDE, TECNICA SIMONA FIORI, TECNICA SIMONA MERLOT, TECNICA SIMONA BIANCO/PLATINO, TECNICA SIMONA GRIGIO/SABBIA, TECNICA SONIA ROSSO, TECNICA SONIA BLU, TECNICA GRAZIA CELESTE, TECNICA GRAZIA SABBIA, TECNICA GRAZIA NERO/GRIGIO, TECNICA CLAUDIA LAMINATO ROSSO, TECNICA KIARA BLU, TECNICA KIARA GHIACCIO PERLA/ARGENTO, TECNICA BEA GRIGIO, TECNICA ALMA BLU, TECNICA CARLA, TECNICA GEA, TECNICA BETTA, TECNICA GAIA LANA ARANCIO, TECNICA GAIA LANA ROSSO, TECNICA GAIA LANA BLU, TECNICA GAIA LANA VERDE, TECNICA ALEXIA.

THE PLANTARE CORK MEMORY CAN BE INSERTED IN: TECNICA CRISTINA ARGENTO, TECNICA ANGELA SABBIA, TECNICA ANGELA NERO, TECNICA MARY BORDEAUX, TECNICA MARY VERDE, TECNICA SARA LILLA, TECNICA SARA FIORI, TECNICA MARY NERO, TECNICA MARY SABBIA, TECNICA MARY BRONZO, TECNICA TANIA, TECNICA LOREN, TECNICA DOROTY NERO, TECNICA DORIS GRIGIO, TECNICA ASIA, TECNICA BENEXA.

THE PLANTARE CORK SLOW CAN BE INSERTED IN: TECNICA RICCARDO NERO, TECNICA RICCARDO, TECNICA RICK NERO, TECNICA RICK, TECNICA BENEXA, TECNICA ANGELO, TECNICA DEVIS, TECNICA DYLAN.

COMFORT
LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

TECNICA NEVA

Certification no. 2710470

Color Ice

Sizes from 35 to 43

Ref. no. EX043600

Diabetic foot risk level

0 1 2 3 see page 4



PADDED
LINING



ADJUSTABLE DOUBLE
VELCRO OPENING



TECNICA VANIA

Certification no. 2710473

Color Gold / Taupe

Sizes from 35 to 43

Ref. no. EX043700

Diabetic foot risk level

0 1 2 3 see page 4



PADDED
LINING



SIDE
INSERT
STRETCH



STRETCH SIDE INSERT
FOR DIFFERENT
PATHOLOGIES



TECNICA MARA VERDE

Certification no. 2755580

Color Green

Sizes from 35 to 43

Ref. no. EX043700

Diabetic foot risk level

0 1 2 3 see page 4



PADDED
LINING



SIDE
INSERT
STRETCH



STRETCH SIDE INSERT





TECNICA ADELE PETROLIO

Certification no. 2421182

Color Petroleum

Sizes from 35 to 43

Ref. no. EX043350

Diabetic foot risk level

① ② ③ see page 4



ADJUSTABLE DOUBLE
VELCRO OPENING



COMFORT LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDLED
LINING

COMFORT
LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

TECNICA ADELE ROSSO

Certification no. 2421181

Color Rosso

Sizes from 35 to 43

Ref. no. EX043350

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA ADELE OCRA

Certification no. 2258636

Color Ocher

Sizes from 35 to 43

Ref. no. EX043350

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA ADELE SUDAN

Certification no. 2258637

Color Sudan

Sizes from 35 to 43

Ref. no. EX043350

Diabetic foot risk level

0 1 2 3 see page 4



COMFORT
LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

TECNICA ADELE FIORI

Certification no. 2421183



Color Flowers

Sizes from 35 to 43

Ref. no. EX043350

Diabetic foot risk level

0 1 2 3 see page 4



ADJUSTABLE DOUBLE
VELCRO OPENING

TECNICA ANNA FIORI

Certification no. 2421184



SIDE
INSERT
STRETCH

Color Flowers

Sizes from 35 to 43

Ref. no. EX043500

Diabetic foot risk level

0 1 2 3 see page 4



STRETCH SIDE INSERT FOR
DIFFERENT PATHOLOGIES

COMFORT
LINEREMOVABLE
INSOLEEXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICEPADDED
LINING**TECNICA SIMONA
FIORI**

Certification no. 2421179

Color Flowers

Sizes from 35 to 43

Ref. no. EX043200

Diabetic foot risk level

0 1 2 3 see page 4

**TECNICA SIMONA
BIANCO/PLATINO**

Certification no. 2710534

Color White/Platinum

Sizes from 35 to 43

Ref. no. EX043100

Diabetic foot risk level

0 1 2 3 see page 4

**TECNICA SONIA ROSSO**

Certification no. 2258640

Color Red

Sizes from 35 to 43

Ref. no. EX043200

Diabetic foot risk level

0 1 2 3 see page 4

**TECNICA SONIA BLU**

Certification no. 2258641

Color Blue

Sizes from 35 to 43

Ref. no. EX043200

Diabetic foot risk level

0 1 2 3 see page 4

**TECNICA SIMONA
GRIGIO/SABBIA**

Certification no. 2710537

Color Grey/Sand

Sizes from 35 to 43

Ref. no. EX043100

Diabetic foot risk level

0 1 2 3 see page 4

**TECNICA GRAZIA
NERO/GRIGIO**

Certification no. 2710526

Color Black/Grey

Sizes from 35 to 43

Ref. no. EX043250

Diabetic foot risk level

0 1 2 3 see page 4



COMFORT
LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

TECNICA CARLOTTA

Certification no. 2755582

Color Green

Sizes from 35 to 43

Ref. no. EX043850

Diabetic foot risk level


0 1 2 3 see page 4



PADDED
LINING



SIDE
INSERT
STRETCH

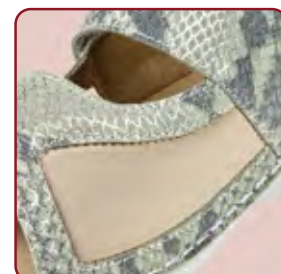

STRETCH SIDE INSERT

TECNICA VITTORIA

Certification no. 2755581


Color Gold

Sizes from 35 to 43

Ref. no. EX043700

Diabetic foot risk level

0 1 2 3 see page 4

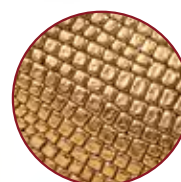

ADJUSTABLE TOTAL
OPENING

TECNICA NINA

Certification no. 2710476

Color Black/Silver

Sizes from 35 to 43

Ref. no. EX044200

Diabetic foot risk level


0 1 2 3 see page 4



PADDED
LINING



SIDE
INSERT
STRETCH

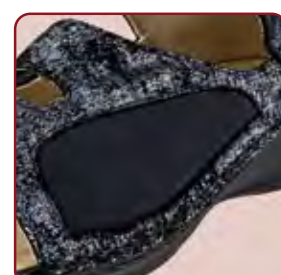

STRETCH SIDE INSERT

COMFORT
LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

TECNICA KIARA BLU

Certification no. 2710515



Color Blue

Sizes from 35 to 43

Ref. no. EX043450

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA KIARA GHIACCIO
PERLA/ARGENTO

Certification no. 2710538



Color Ice/Pearl/
Silver

Sizes from 35 to 43

Ref. no. EX043450

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA BEA GRIGIO

Certification no. 2258643



Color Grey

Sizes from 35 to 43

Ref. no. EX043550

Diabetic foot risk level

0 1 2 3 see page 4

COMFORT
LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

TECNICA CRISTINA
ARGENTO

Certification no. 2258639



Color Silver

Sizes from 35 to 42

Ref. no. EX043850

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA ANGELA SABBIA

Certification no. 2710555



Color Sand

Sizes from 35 to 42

Ref. no. EX043850

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA ANGELA NERO

Certification no. 2710525



Color Black

Sizes from 35 to 42

Ref. no. EX043850

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA MARY
BORDEAUX

Certification no. 2421191



Color Bordeaux

Sizes from 35 to 42

Ref. no. EX043700

Diabetic foot risk level

0 1 2 3 see page 4

COMFORT
LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING



TECNICA ZARA

Certification no. 2710478

Color White / Copper

Sizes from 35 to 43

Ref. no. EX043850

Diabetic foot risk level

0 1 2 3 see page 4



PADDED
LINING



SIDE
INSERT
STRETCH



TECNICA EDEN

Certification no. 2710477

Color Luxor

Sizes from 35 to 43

Ref. no. EX043850

Diabetic foot risk level

0 1 2 3 see page 4



PADDED
LINING



SIDE
INSERT
STRETCH



COMFORT
LINEREMOVABLE
INSOLEEXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICEPADDED
LINING**TECNICA GEMMA**

Certification no. 2710483

Color Sand

Sizes from 35 to 42

Ref. no. EX043800

Diabetic foot risk level

0 1 2 3 see page 4

PADDED
LININGSIDE
INSERT
STRETCH**TECNICA FLORA**

Certification no. 2755583

Color Green

Sizes from 35 to 42

Ref. no. EX043700

Diabetic foot risk level

0 1 2 3 see page 4

PADDED
LININGSIDE
INSERT
STRETCH

WOMAN

100% MADE IN ITALY

COMFORT LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

TECNICA MARY NERO

Certification no. 2710553



SIDE
INSERT
STRETCH

Color Black

Sizes from 35 to 42

Ref. no. EX043700

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA MARY SABBIA

Certification no. 2710550



SIDE
INSERT
STRETCH

Color Sand

Sizes from 35 to 42

Ref. no. EX043700

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA MARY BRONZO

Certification no. 2710509



SIDE
INSERT
STRETCH

Color Bronze

Sizes from 35 to 42

Ref. no. EX043700

Diabetic foot risk level

0 1 2 3 see page 4



**TECNICA ALMA BLU**

Certification no. 2710564

Color Blue

Sizes from 35 to 42

Ref. no. EX043450

Diabetic foot risk level

0 ① ② ③ see page 4



REAL WOOL INSIDE

**COMFORT
LINE**REMOVABLE
INSOLEEXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICEPADDED
LINING

COMFORT LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

TECNICA CARLA

Certification no. 2421185

Color Grey

Sizes from 35 to 43

Ref. no. EX043900

Diabetic foot risk level

① ② ③ see page 4



STRETCH
UPPER



ADJUSTABLE VELCRO
OPENING

TECNICA GEA

Certification no. 2421186

Color Beige

Sizes from 35 to 43


Ref. no. EX043900

Diabetic foot risk level

① ② ③ see page 4



STRETCH
UPPER

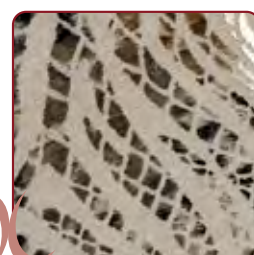

TECNICA BETTA

Certification no. 2421187

Color Blue

Sizes from 35 to 43


Ref. no. EX043900

Diabetic foot risk level

① ② ③ see page 4



STRETCH
UPPER

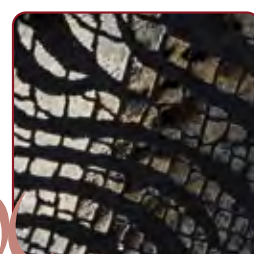

COMFORT
LINEREMOVABLE
INSOLEEXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICEPADDED
LINING

TECNICA LOREN

Certification no. 2421197



Color Black

Sizes from 35 to 42

Ref. no. EX043650

Diabetic foot risk level

0 ① ② ③ see page 4

STRETCH
UPPERTECNICA DOROTY
NERO

Certification no. 2710567



Color Black

Sizes from 35 to 42

Ref. no. EX043650

Diabetic foot risk level

0 ① ② ③ see page 4

STRETCH
UPPER

TECNICA ASIA

Certification no. 2710507

Color Black

Sizes from 35 to 42

Ref. no. EX043650

Diabetic foot risk level

0 ① ② ③ see page 4

STRETCH
UPPERINSOLE AND UPPER
LINING COVERED WITH ON
STEAM® MICROFIBER

COMFORT
LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

**TECNICA RICCARDO
NERO**

Certification no. 2418173

Color Black

Sizes from 39 to 50

Ref. no. EX043950

Diabetic foot risk level
0 1 2 3 see page 4



TECNICA RICCARDO

Certification no. 2710556

Color Brown

Sizes from 39 to 50

Ref. no. EX043950

Diabetic foot risk level
0 1 2 3 see page 4



TECNICA RICK NERO

Certification no. 2418177

Color Black

Sizes from 39 to 50

Ref. no. EX044300

Diabetic foot risk level
0 1 2 3 see page 4



TECNICA RICK

Certification no. 2710504

Color Brown

Sizes from 39 to 50

Ref. no. EX044300

Diabetic foot risk level
0 1 2 3 see page 4



TECNICA BENEXA

Certification no. 2710508

Color Black

Sizes from 35 to 50

Ref. no. EX044650

Diabetic foot risk level
0 1 2 3 see page 4



UNISEX

COMFORT
LINE



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



PADDED
LINING

TECNICA DEVIS

Certification no. 2421194



Color Blue

Sizes from 35 to 50

Ref. no. EX042950

Diabetic foot risk level

0 1 2 3 see page 4



ADJUSTABLE VELCRO
CLOSURE

TECNICA DYLAN

Certification no. 2421195




Color Beige

Sizes from 35 to 50

Ref. no. EX042950

Diabetic foot risk level

0 1 2 3 see page 4

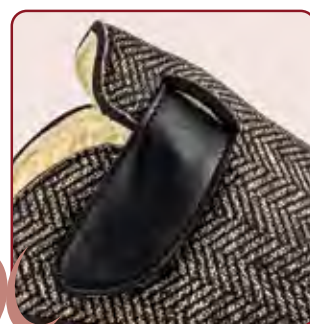

REAL WOOL INSIDE
AND ADJUSTABLE
VELCRO CLOSURE



REHABILITATION LINE



REHABILITATION LINE

COMFORT, PROTECTION, AND VERSATILITY FOR EVERYDAY FOOT WELL-BEING

The Rehabilitation Line is designed to offer lightweight, comfortable, and technically advanced footwear, specifically intended for home use and for the needs of patients in the post-operative phase, elderly individuals, or people with chronic foot conditions. Characterized by an extra-wide fit (fitting 14) and an adjustable structure with Velcro closures, this line is ideal for managing diabetic foot. However, it also stands out for its multifunctional versatility, making it effective in the daily care of many other clinical conditions.

Indeed, the Rehabilitation Line represents a valid option for patients suffering from: hallux valgus, hammer toes, metatarsalgia, rheumatic foot, edema, post-traumatic or post-surgical deformities, and painful conditions affecting the plantar area and/or the forefoot. Each model is designed to ensure maximum facility of wear, even in the presence of bandages or swelling, and can accommodate custom orthotics up to 9 mm—providing a functional base for personalized treatments.

TECHNICAL FOCUS: OUTSOLE MODULARITY BASED ON THE CLINICAL RISK OF DIABETIC FOOT

A distinctive feature of this line is the ability to customize the outsole according to the patient's clinical profile and the rehabilitation phase:

- **FLEXIBLE OUTSOLE**, for safe and comfortable use in the home environment and in the absence of ulcers or structural issues.
- **SEMI-RIGID OUTSOLE**, designed to improve biomechanics and support the foot in moderate-risk conditions.
- **RIGID OUTSOLE**, for the treatment of high-risk diabetic foot and for advanced offloading and containment needs.

This modularity makes the footwear highly adaptable, allowing healthcare professionals to customize the therapeutic approach based on the pathology and the environment of use.

Although designed for indoor use, the Rehabilitation Line offers a complete technical structure, featuring supportive heel counters, self-molding fabrics, and a standard breathable, antibacterial S-Air insole.

It is, in every respect, a high-performance medical shoe designed to ensure lightweight comfort, breathability, and facility of use throughout the day.

The Rehabilitation Line is a modern, customizable, and versatile therapeutic solution, capable of providing maximum support in the home setting to patients with diverse clinical needs—without compromising comfort, functionality, or protection.

THE DAILY CHOICE FOR REHABILITATION AND PREVENTION

PLANTARE S-AIR FISIO

Certification no. 2702558

Sizes from 35 to 50

Ref. no. EX010780



PLANTARE S-RIGID FISIO

Certification no. 2702559

Sizes from 35 to 50

Ref. no. EX011710



PLANTARE S-CLOUD FISIO

Certification no. 2142383

Sizes from 35 to 50

Ref. no. EX012050



PLANTARE S-SEMIRIGID FISIO


Certification no. 2798981

Sizes from 35 to 50

Ref. no. EX011710



REHABILITATION LINE FIT 11

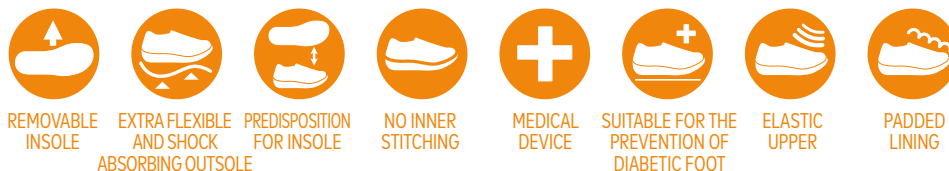

TECNICA MARISA

Certification no. 2470727



Color Cream

Sizes from 35 to 42

Ref. no. EX043950

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA MARTINA

Certification no. 2470728


Color Cream

Sizes from 35 to 42

Ref. no. EX043950

Diabetic foot risk level

0 1 2 3 see page 4

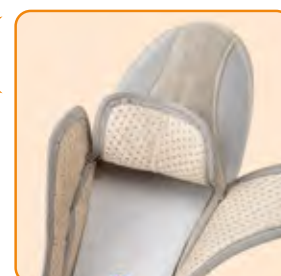

VELCRO OPENING

TECNICA MIRIAM

Certification no. 2470729



Color Black


Sizes from 35 to 42

Ref. no. EX043950

Diabetic foot risk level

0 1 2 3 see page 4

REHABILITATION LINE FIT 14

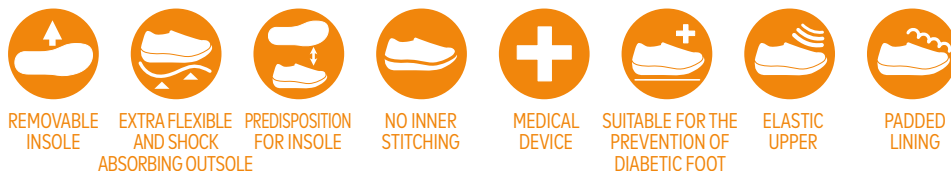

TECNICA AGRI

Certification no. 2470730




Color Night blue

Sizes from 35 to 50

Ref. no. EX043950

Diabetic foot risk level

0 1 2 3 see page 4

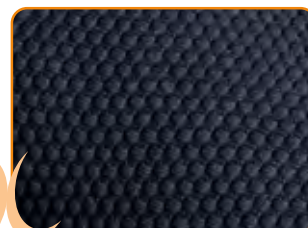

TOMAIA TRASPIRANTE

TECNICA ANDRIA

Certification no. 2470731




Color Green

Sizes from 35 to 50

Ref. no. EX043950


Diabetic foot risk level

0 1 2 3 see page 4

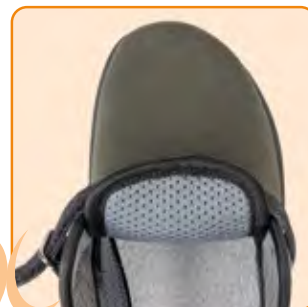

VELCRO OPENING
REMOVABLE FOOTBED

REHABILITATION LINE FIT 14

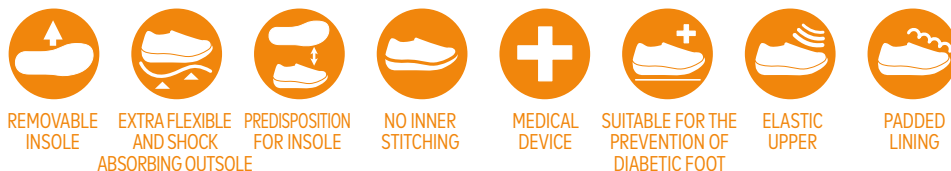

TECNICA AVOLA

Certification no. 2702577



Color	Black
Sizes	from 35 to 50
Ref. no.	EX045450
Diabetic foot risk level	
0 1 2 3 see page 4	

TECNICA DYNAMIC UOMO

Certification no. 2702574



Color	Black
Sizes	from 39 to 50
Ref. no.	EX013550
Diabetic foot risk level	
0 1 2 3 see page 4	

UOMO


TECNICA GOLD 1

Certification no. 2702538



Color	Black
Sizes	from 35 to 50
Ref. no.	EX013900
Diabetic foot risk level	
0 1 2 3 see page 4	

REHABILITATION LINE FIT 14

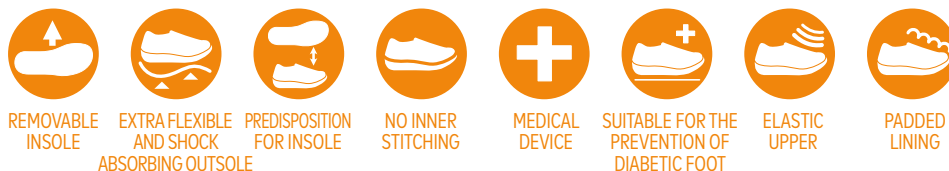

TECNICA FLEX

Certification no. 2702563



Color Grey/Black

Sizes from 35 to 50

Ref. no. EX013350

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA DYNAMIC AIR

Certification no. 2702570



Color Black

Sizes from 35 to 50

Ref. no. EX012550

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA SOFT

Certification no. 2702576


Color Black

Sizes from 35 to 50


Ref. no. EX012550

Diabetic foot risk level

0 1 2 3 see page 4

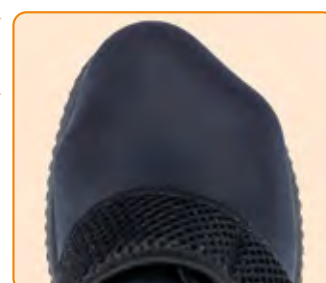

REHABILITATION LINE FIT 14

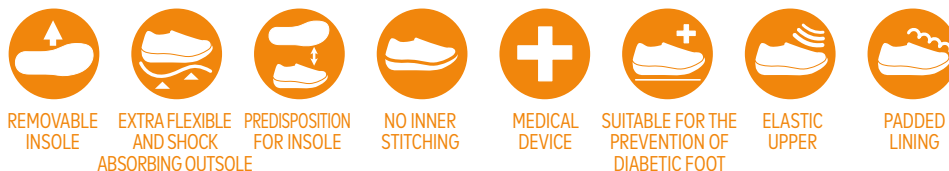

TECNICA 3T

Certification no. 2702543



Color Grey

Sizes from 35 to 50

Ref. no. EX013000

Diabetic foot risk level

0 1 2 3 see page 4



CHIUSURA VELCRO E
APERTURA TOTALE

TECNICA 3S

Certification no. 2702553




Color Black

Sizes from 35 to 50


Ref. no. EX013200

Diabetic foot risk level

0 1 2 3 see page 4

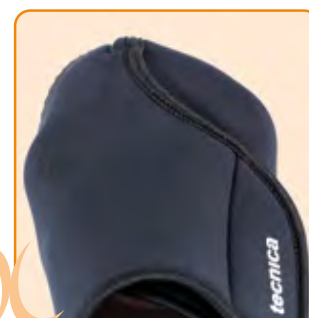

REHABILITATION
LINE
FIT 14

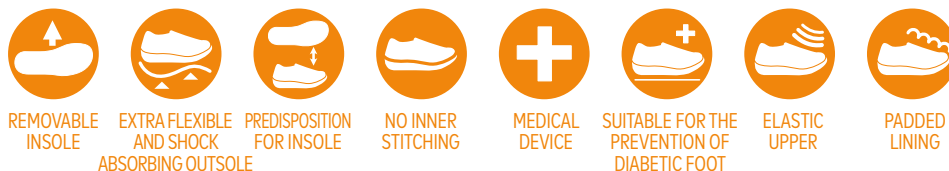

TECNICA 3M

Certification no. 2702544



Color Blue

Sizes from 35 to 50

Ref. no. EX013300

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA 3E

Certification no. 2702545



Color Black

Sizes from 35 to 50

Ref. no. EX013200

Diabetic foot risk level

0 1 2 3 see page 4

TECNICA 3A

Certification no. 2702542

Color Grey

Sizes from 35 to 50

Ref. no. EX013200

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA 1

Certification no. 2702539

Color Grey

Sizes from 35 to 50


Ref. no. EX013500

Diabetic foot risk level

0 1 2 3 see page 4



REHABILITATION LINE FIT 14

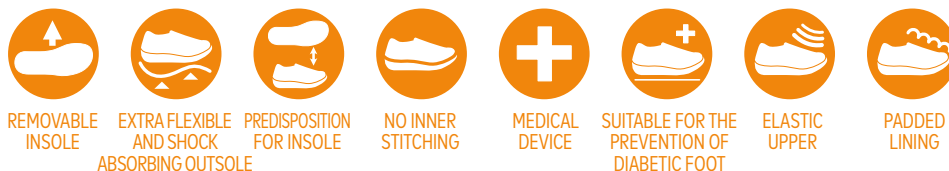

TECNICA 3 KLE BLU

Certification no. 2702573

Color Blue

Sizes from 35 to 50

Ref. no. EX013350

Diabetic foot risk level

0 1 2 3 see page 4



Trattamento antibatterico agli ioni d'argento Silver Plus resistente all'acqua e non nocivo sulla pelle.

TECNICA 3 AGNELLO

Certification no. 2702555

Color Blue

Sizes from 35 to 50

Ref. no. EX013300

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA 3 A KLE NERO

Certification no. 2702572

Color Black

Sizes from 35 to 50

Ref. no. EX013450

Diabetic foot risk level

0 1 2 3 see page 4



Trattamento antibatterico agli ioni d'argento Silver Plus resistente all'acqua e non nocivo sulla pelle.

REHABILITATION LINE FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



ELASTIC
UPPER



PADDED
LINING

TECNICA 4T

Certification no. 2702546

Color	Grey
Sizes	from 35 to 50
Ref. no.	EX012500

Diabetic foot risk level
0 1 2 3 see page 4



TECNICA 4M

Certification no. 2702547

Color	Blue
Sizes	from 35 to 50
Ref. no.	EX012500

Diabetic foot risk level
0 1 2 3 see page 4



TECNICA 11

Certification no. 2702550

Color	Grey
Sizes	from 35 to 50
Ref. no.	EX012750

Diabetic foot risk level
0 1 2 3 see page 4



TECNICA 10

Certification no. 2702568

Color	Grey
Sizes	from 35 to 50
Ref. no.	EX012150

Diabetic foot risk level
0 1 2 3 see page 4



TECNICA T2 AIR


Certification no. 2702569

Color	Grey
Sizes	from 35 to 50
Ref. no.	EX013000

Diabetic foot risk level
0 1 2 3 see page 4



REHABILITATION LINE FIT 14

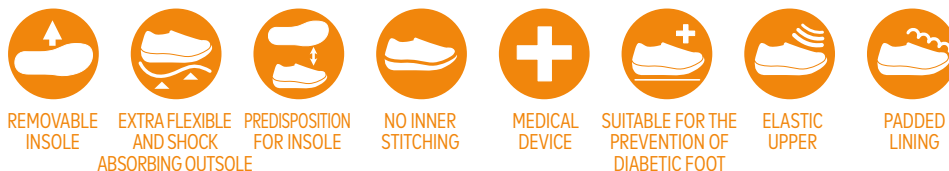

TECNICA 5T

Certification no. 2702564

Color Grey

Sizes from 35 to 50

Ref. no. EX013000

Diabetic foot risk level

0 1 2 3 see page 4



TOTAL OPENING

TECNICA 7

Certification no. 2702567

Color Grey

Sizes from 35 to 50

Ref. no. EX012850

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA 6

Certification no. 2702548

Color Grey

Sizes from 35 to 50

Ref. no. EX012350

Diabetic foot risk level

0 1 2 3 see page 4



BAMBOO REHABILITATION LINE

MADE OF AN UPPER COVERED IN BAMBOO, A GIFT OF NATURE FOR FOOT WELL-BEING. BAMBOO IS NATURAL AND GENTLE ON THE SKIN, PREVENTS IRRITATIONS AND BAD ODORS, AND IS SUSTAINABLE AND ENVIRONMENTALLY FRIENDLY!



NATURAL ANTIBACTERIAL

Bamboo contains a natural antimicrobial agent called bamboo kun, which inhibits the growth of bacteria and fungi. This property makes bamboo ideal for direct skin contact applications, helping to prevent irritations and bad odors.



MOISTURE ABSORPTION

Bamboo has an extraordinary ability to absorb moisture, up to 3-4 times more than cotton. This helps to keep the skin dry, reducing the risk of rubbing or irritation.



BREATHABILITY AND THERMOREGULATION

Bamboo fibers are highly breathable and allow a natural temperature regulation. They keep the feet cool in summer and warm in winter, ensuring optimal comfort throughout the year.



NATURAL SOFTNESS

Bamboo fibers are incredibly soft, similar to silk, but much more durable.



HYPOALLERGENIC

Thanks to its natural composition and the absence of harsh chemicals, bamboo is an excellent choice for those with allergies or sensitive skin.




BAMBOO IS ONE OF THE MOST ENVIRONMENTALLY FRIENDLY NATURAL RESOURCES IN THE WORLD: IT GROWS RAPIDLY WITHOUT THE NEED FOR PESTICIDES OR FERTILIZERS. IT REGENERATES NATURALLY AFTER CUTTING. IT ABSORBS LARGE AMOUNTS OF CARBON DIOXIDE AND PRODUCES OXYGEN.



UPPER COVERED WITH BAMBOO FIBER



REHABILITATION LINE FIT 14

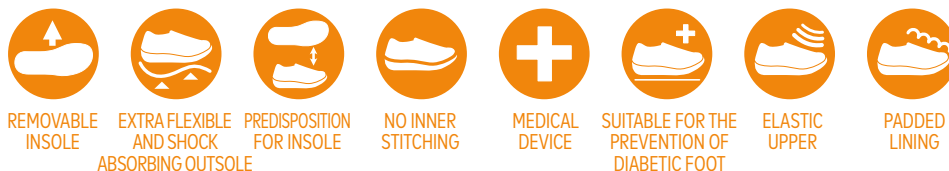

TECNICA 3 JAK GRIGIO

Certification no. 2755197

Color Grey

Sizes from 35 to 50

Ref. no. EX013100

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA 3L JAK BLU

Certification no. 2755198


Color Blue

Sizes from 35 to 50

Ref. no. EX013100

Diabetic foot risk level

0 1 2 3 see page 4

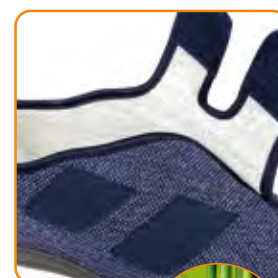

TECNICA 3A SPIG NERO

Certification no. 2755199

Color Black

Sizes from 35 to 50

Ref. no. EX013300

Diabetic foot risk level

0 1 2 3 see page 4



REHABILITATION LINE FIT 14



REMOVABLE
INSOLE



EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



ELASTIC
UPPER



PADDED
LINING

TECNICA 11 JAK GRIGIO

Certification no. 2755200

Color Grey

Sizes from 35 to 50

Ref. no. EX012850

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA 11 JAK BEIGE

Certification no. 2755201

Color Beige

Sizes from 35 to 50

Ref. no. EX012850

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA 11 SPIG NERO

Certification no. 2755202

Color Black

Sizes from 35 to 50


Ref. no. EX012850

Diabetic foot risk level

0 1 2 3 see page 4



REHABILITATION LINE FIT 14

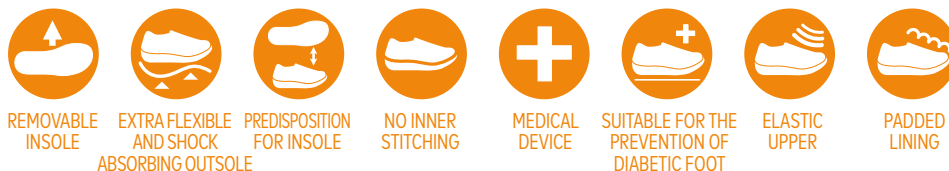

REMOVABLE
INSOLE

EXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLE

PREDISPOSITION
FOR INSOLE

NO INNER
STITCHING

MEDICAL
DEVICE

SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT

ELASTIC
UPPER

PADDED
LINING

TECNICA 10 JAK GRIGIO

Certification no. 2755203

Color Grey

Sizes from 35 to 50

Ref. no. EX012250

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA 10 JAK BEIGE

Certification no. 2755204

Color Beige

Sizes from 35 to 50

Ref. no. EX012250

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA 6 JAK BLU

Certification no. 2755205

Diabetic foot risk level

0 1 2 3 see page 4

Color Blue

Sizes from 35 to 50

Ref. no. EX012450



TECNICA 4 JAK GRIGIO

Certification no. 2755206

Diabetic foot risk level

0 1 2 3 see page 4


Color Grey

Sizes from 35 to 50

Ref. no. EX012600



REHABILITATION LINE XXL VOLUME

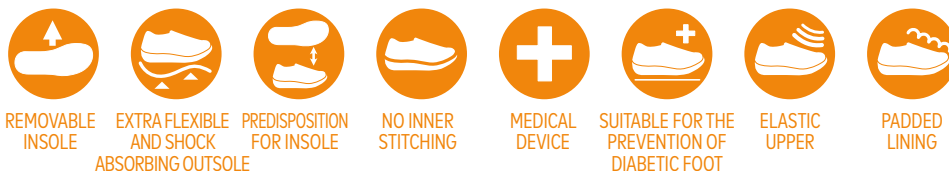REMOVABLE
INSOLEEXTRA FLEXIBLE
AND SHOCK
ABSORBING OUTSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICESUITABLE FOR THE
PREVENTION OF
DIABETIC FOOTELASTIC
UPPERPADDED
LINING

XXL REHABILITATION FOOTWEAR

THIS NEW LINE OF XXL REHABILITATION FOOTWEAR WITH EXTRA-LARGE VOLUME IS DESIGNED TO PROVIDE MAXIMUM COMFORT AND SUPPORT FOR THOSE WHO NEED SHOES WITH GREAT SPACE.

Ideal for the treatment of foot pathologies, deformities, rheumatic foot, metatarsalgia, diabetic foot (risk levels 0 and 1), and post-operative rehabilitation.

MAIN FEATURES

XXL VOLUME

Designed to provide extra space and comfort, even in the presence of deformities or swelling.

ADJUSTABLE VELCRO CLOSURE

Allows a perfect and customizable fitting.

DESIGNED TO ACCOMMODATE 8 MM INSOLES

Enables the insertion of custom orthotic insoles.

COTTON LINING TREATED WITH SILVER IONS

Provides effective antibacterial treatment, helping to keep the foot dry and fresh.

SLIP-RESISTANT OUTSOLE

Ensures stability and safety during walking.

This footwear is the ideal solution for those seeking a perfect balance of comfort, protection, and functionality during the rehabilitation process or when dealing with specific foot needs.

TECNICA 13 XXL NERO

Certification no. 2755208

Color Black

Sizes from 36 to 50

Ref. no. EX013900

Diabetic foot risk level

0 1 2 3 see page 4



TECNICA 11 XXL NERO

Certification no. 2755207

Color Black

Sizes from 36 to 50


Ref. no. EX013500

Diabetic foot risk level


0 1 2 3 see page 4



REHABILITATION LINE

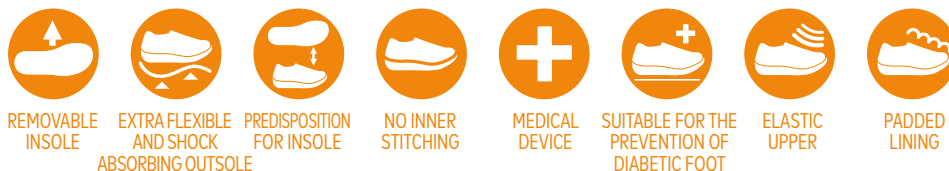

INTERNAL MEASURES

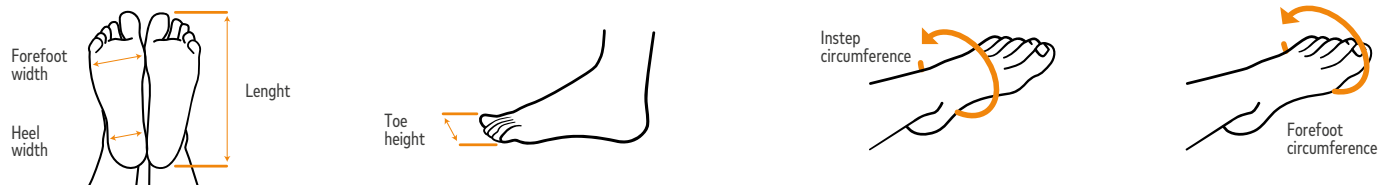

MEASURES OF THE MODELS: AGRI, ANDRIA, AVOLA, DYNAMIC UOMO, GOLD 1, FLEX, DYNAMIC AIR, SOFT, 3T, 3S, 3M, 3E, 3A, 1, 3 KLE BLU, 3 AGNELLO, 3A KLE NERO, 4T, 4M, 11, 10, T2 AIR, 5T, 7, 6, 3 JAK GRIGIO, 3L JAK BLU, 3A SPIG NERO, 11 JAK GRIGIO, 11 JAK BEIGE, 11 SPIG NERO, 10 JAK GRIGIO, 10 JAK BEIGE, 6 JAK BLU, 4 JAK GRIGIO.

FIT 14															
Sizes	UM	35	36	37	38	39	40	41	42	43	44	45	46	47/48	49/50
Lenght	mm	225	232	239	245	252	259	265	272	279	285	292	299	308	323
Forefoot width	mm	83	84,5	86	87,5	89	90,5	92	93,5	95	96,5	98	99,5	102	105
Heel width	mm	63,5	65	66	67	68	69,5	71	72	73	74	75	77	78,5	81
Forefoot circumference	mm	235	239,5	244	249	253,5	258	262,5	267	271,5	276	280,5	285	292	301
Toe height	mm	35	36	36,5	37	37	37,5	38,5	39	39,5	40	41	41,5	42,5	43,5
Instep circumference	mm	346	357	366	374	383	387	396	403	413	422	431	439	451	470

MEASURES OF THE MODELS: 13 XXL NERO, 11 XXL NERO.

XXL REHABILITATION MODEL FIT														
Sizes	UM	36	37	38	39	40	41	42	43	44	45	46	47/48	49/50
Lenght	mm	235	242	249	255	262	269	275	282	289	295	302	315	329
Forefoot width	mm	99	100,5	102	103,5	105	106,5	108	109,5	111	112,5	114	117	120
Heel width	mm	69	70,5	72	73,5	75	76,5	78	79,5	81	82,5	84	87	90
Forefoot circumference	mm	281	286	290	294	298	303	307	311	315	320	324	332	340
Toe height	mm	35	35,5	36	36,5	37	37,5	38	38,5	39	39,5	40	41	42



DIABETIC FOOT BRACES AND POST-OPERATIVE DEVICES LINE



DIABETIC FOOT BRACES AND POST-OPERATIVE DEVICES LINE

ADVANCED SOLUTIONS FOR THE TREATMENT OF DIABETIC FOOT AND POST-OPERATIVE RECOVERY.

The line dedicated to foot-ankle braces and post-operative care is designed to meet the most complex clinical needs, particularly in cases of acute diabetic foot, neuropathic and

ischemic ulcers of the forefoot and midfoot, and for post-surgical treatment following foot surgeries.

These technical shoes are designed to provide offloading, protection, and stability throughout all phases of treatment and recovery, both at home and in clinical settings.

MAIN TECHNICAL FEATURES:

→ Super-soft padded upper, seamless inside:

The interior of the shoe is designed to be as gentle as possible. Thanks to soft materials and the absence of seams, this solution reduces the risk of rubbing, abrasions, or trauma, making it ideal even in the presence of lesions, fragile skin, or bandages.

→ Multi-velcro closure for adjustable volume:

Multiple adjustable straps allow the shoe to perfectly adapt to the foot, even in cases of swelling, bandages, or structural deformities. The full opening also facilitates putting on the shoe, even for patients with limited mobility.

→ Thermoformable and stabilizing heel counter:

The rigid heel counter can be thermoformed to provide customized containment, ensuring effective stabilization of the ankle and rearfoot during walking.

→ Rocking rigid outsole with slip-resistant insert:

Designed to promote a rolling gait, it helps reduce pressure on the forefoot, protecting areas at risk of ulceration. The central slip-resistant insert ensures excellent grip and safety while walking.

→ Lightweight and balanced structure:

Despite its advanced technical features, the models maintain a lightweight structure and functional design, making daily use easier without causing patient fatigue.



THE FOOT BRACES AND POST-OPERATIVE LINE REPRESENTS THE PERFECT COMBINATION OF TECHNOLOGY, THERAPEUTIC FUNCTIONALITY, AND PROTECTIVE COMFORT, OFFERING HEALTHCARE PROFESSIONALS AN EFFECTIVE AND SAFE SOLUTION TO SUPPORT PATIENTS THROUGH THE MOST DELICATE PHASES OF THEIR TREATMENT PROCESS.

MAIN APPLICATIONS:

- Treatment and prevention of plantar ulcers in diabetic patients
- Post-operative offloading of the foot and forefoot
- Post-surgical management of toes, metatarsals, and midfoot
- Presence of bandages, edema, or deformities
- Need for selective offloading and containment of the foot and ankle



Thermoformable and stabilizing heel counter

Super-soft padded upper with seamless interior

Lightweight and balanced structure

Multi-velcro closure for adjustable volume

Rigid Rocker outsole with slip-resistant insert



DIABETIC FOOT BRACES AND POST-OPERATIVE DEVICES LINE

VARIABLE VOLUME - AMBIDEXTROUS

REMOVABLE
INSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICESUITABLE FOR THE
PREVENTION OF
DIABETIC FOOTELASTIC
UPPERPADDED
LININGROCKING
OUTSOLE

HELIOS REHA

Certification no. 2702880



Color Black

Sizes from S to XXL

Ref. no. EX027900



REMOVABLE REAR COVER

INDICATIONS

The rigid brace is indicated for immobilization and rehabilitation of the foot after surgery:

- Achilles tendon surgery

CHARACTERISTICS

The rigid brace is a medical device that is adjustable using hook-and-loop straps, ambidextrous, and made with two lateral support bars, as well as a in height-adjustable heel. The back part is open to prevent pressure points in sensitive areas.

Helios Reha is used to:

- Prevent ankle movement and the resulting extension of the Achilles tendon.
- Allow for programmed rehabilitation after Achilles tendon surgery (adjustable heel).
- Avoid pressure points in sensitive areas through its rear opening.

Sizes				
S	M	L	XL	XXL
34-35-36	37-38-39	40-41-42	43-44-45	46-47-48

HELIOS HEEL

Certification no. 2702881



Color Black

Sizes from S to XXL

Ref. no. EX027700



REMOVABLE REAR COVER

INDICATIONS

- Management of heel ulcers (both plantar and dorsal)
- Achilles tendon ulcer/wound
- Heel decubitus ulcers (both plantar and dorsal)
- Post-operative rehabilitation of the heel

CHARACTERISTICS

The HELIOS HEEL brace is a medical device that is adjustable using hook-and-loop straps, ambidextrous, and composed of two lateral reinforcement bars, as well as a in height-adjustable heel. The back is open to avoid pressure points in sensitive areas.

Sizes				
S	M	L	XL	XXL
34-35-36	37-38-39	40-41-42	43-44-45	46-47-48

DIABETIC FOOT BRACES AND POST-OPERATIVE DEVICES LINE

VARIABLE VOLUME - AMBIDEXTROUS

REMOVABLE
INSOLEPREDISPOSITION
FOR INSOLENO INNER
STITCHINGMEDICAL
DEVICESUITABLE FOR THE
PREVENTION OF
DIABETIC FOOTELASTIC
UPPERPADDED
LININGROCKING
OUTSOLE

TECNICA DIAB

Certification no. 2702877



Color Black

Sizes from S to XXL

Ref. no. EX028200



TECNICA CT-DIAB

INDICATIONS

- Diabetic plantar ulcer
- Plantar ulcer
- Post-amputation rehabilitation (up to trans-metatarsal amputation)
- Post-operative rehabilitation of the tarsus
- Fracture of the tarsus
- Management of the chronic phase of Charcot's foot
- Malleolus fracture

CHARACTERISTICS

High leg-foot brace, which is adaptable to the patient, allows the displacement of part of the load on the ankles and calves. The outsole, thanks to the features of movement control and pressure redistribution decompresses the areas affected by the injury, improving thus significantly the ambulation. High thermoformable stiffener. Rigid, anti-flexion and rocking outsole. Variable volume.

Predisposition for the S-Cloud Diab insole.

Predisposition for the S-Tal Diab insole.

Sizes				
S	M	L	XL	XXL
34-35-36	37-38-39	40-41-42	43-44-45	46-47-48

TECNICA POST

Certification no. 2702885



Color Black

Sizes from S to XXL

Ref. no. EX026200



TECNICA CT-DIAB

INDICATIONS

- Diabetic plantar ulcer with vascular complications affecting the leg (combined with the S-Cloud Diab insole)
- Neuroischemic ulcer
- Management of the injured foot or foot with ulcerative risk
- Compound fracture of the metatarsals
- Post-amputation rehabilitation of the toes
- Fracture of the toes
- Surgical treatment of the toes
- Pre-surgery stabilization

CHARACTERISTICS

Low leg-foot brace, which is adaptable to the patient. The outsole, thanks to the features of movement control and pressure redistribution decompresses the areas affected by the injury, improving thus significantly the ambulation and the recovery time.

High thermoformable stiffener.

Rigid, anti-flexion and rocking outsole. Variable volume.

Predisposition for the S-Cloud Diab insole.

Predisposition for the S-Tal Diab insole.

Sizes				
S	M	L	XL	XXL
34-35-36	37-38-39	40-41-42	43-44-45	46-47-48

DIABETIC FOOT BRACES AND POST-OPERATIVE DEVICES LINE

VARIABLE VOLUME - AMBIDEXTROUS



REMOVABLE
INSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



ELASTIC
UPPER



PADDED
LINING



ROCKING
OUTSOLE

TECNICA TERAP

Certification no. 2702886



Color Black

Sizes from S to XXL

Ref. no. EX025900

INDICATIONS

- Diabetic ulcerations
- Plantar ulcerations
- Neuroischemic ulcerations
- Fracture of the metatarsals
- Surgery and fractures of the toes
- Rehabilitation after toes amputation
- Pre-surgery stabilization
- Instability of the hindfoot

CHARACTERISTICS

Easy fit and maximum safety for the patient during the acute phase. The thermoformable upper in pel-flex, which is partially waterproof, guarantees thermal and hygienic protection to the foot. Radiolucent and anti-friction. Thermoformable stiffener.

Rigid, anti-flexion and rocking outsole. Variable volume. Predisposition for the insertion of the S-Cloud Diab insole.

Sizes				
S	M	L	XL	XXL
34-35-36	37-38-39	40-41-42	43-44-45	46-47-48

TECNICA HEEL

Certification no. 2702896



Color Black

Sizes from S to XXL

Ref. no. EX023400

INDICATIONS

- Post-surgery rehabilitation of the heel area
- Wounds or ulcers of the heel area
- Post-surgery fractures and/or calcaneus ulcers (second treatment) and ulcers at the insertion of the Achilles tendon
- Bursitis

CHARACTERISTICS

The balance of the weight is moved to the front of the foot thanks to the unloading insole and thanks to the outsole, which is designed to totally release the heel from the pressure load. The square shape of the outsole offers more space for bandages. Movable rear-cover to preserve the heel from impacts. Non-slip outsole to ensure maximum safety during walking. Total opening with velcro-straps. Removable insole.

Sizes			
S	M	L	XL
34-35-36	37-38-39	40-41-42	43-44-45

TECNICA CT DIAB

Certification no. 2702893

Color Black

Sizes from S to XXL

Ref. no. EX021370



PLANTARE S-CLOUD DIAB

Certification no. 2702901

Sizes from S to XXL

Ref. no. EX021420



PLANTARE S-TAL DIAB

Certification no. 2702895

Sizes from S to XXL

Ref. no. EX021250



DIABETIC FOOT BRACES AND POST-OPERATIVE DEVICES LINE

VARIABLE VOLUME - AMBIDEXTROUS



REMOVABLE
INSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



ELASTIC
UPPER



PADDED
LINING



ROCKING
OUTSOLE

TECNICA JOLLY

Certification no. 2702887



Color Black

Sizes from S to XXL

Ref. no. EX023900



TECNICA CT-DIAB

INDICATIONS

- Post-surgery of hallux valgus
- Hammer toes
- Morton's neuroma
- Fracture of the toes
- Nail removal
- Plantar wart removal

CHARACTERISTICS

Rehabilitative foot orthosis, with radiolucent upper, which is antifriction and washable. The non-slip, rigid and rocking outsole unloads the pressure from the forefoot, giving thus greater stability without altering the posture, ensuring a rapid recovery of all functions. Thermoformable stiffener. Variable volume. Predisposition for the S-Cloud Diab insole. Predisposition for the S-Tal Diab insole.

Sizes				
S	M	L	XL	XXL
34-35-36	37-38-39	40-41-42	43-44-45	46-47-48

TECNICA VALGUS

Certification no. 2702888



Color Black

Sizes from S to XXL

Ref. no. EX023200

INDICATIONS



- Post-surgery of hallux valgus, hammer toes, Neuroma of Morton
- Fracture of the toes, nail removal
- Plantar wart removal
- Best solution for diabetic feet and indicated in presence of trans-metatarsal amputations
- Ray amputations
- Chronic Charcot of the forefoot and of the midfoot
- Minor amputations of the forefoot

CHARACTERISTICS

The rigid and rocking outsole cushions the step and immobilizes the metatarsophalangeal area, reducing the pressure, allowing thus a total discharge of the trans-metatarsal area. Allows to accommodate deformities and allows to avoid dangerous frictions. The wide plant and the anti-stumbling outsole ensure maximum stability during walking. Higher stabilizing stiffener for better foot support. Variable volume. Removable insole.

Sizes				
S	M	L	XL	XXL
34-35-36	37-38-39	40-41-42	43-44-45	46-47-48

INTERNAL MEASURES

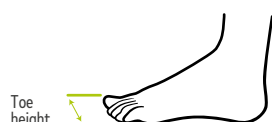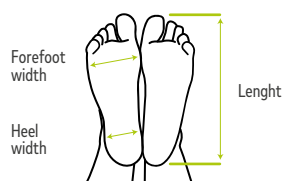

MEASURES OF THE MODELS: DIAB, POST, TERAP, JOLLY, VALGUS.

Sizes	UM	S	M	L	XL	XXL
Length	mm	230	250	270	290	310
Forefoot width	mm	95	100	105	110	115
Heel width	mm	78	80	82	86	90

DIABETIC FOOT BRACES AND POST-OPERATIVE DEVICES LINE

VARIABLE VOLUME - AMBIDEXTROUS



REMOVABLE
INSOLE



PREDISPOSITION
FOR INSOLE



NO INNER
STITCHING



MEDICAL
DEVICE



SUITABLE FOR THE
PREVENTION OF
DIABETIC FOOT



ELASTIC
UPPER



PADDED
LINING



ROCKING
OUTSOLE

TECNICA STELT

Certification no. 2702902

Color Black/White

Sizes from 35 to 48

Ref. no. EX023700



INDICATIONS

- Rigid biomechanical outsole
- 20mm outsole which minimizes discrepancy with the contralateral limb
- S-Tal Stelt insole included
- Stabilizing and thermoformable stiffener
- Multi-velcro closure
- Ambidextrous
- Predisposition for the insertion of the S-Cloud Stelt unloading insole

CHARACTERISTICS

Tecnica STELT is an ambidextrous footwear, indicated both as a brace and as a post-operative device. It is particularly suitable for the relieve of pressure in case of neuropathic and ischemic ulcers in the forefoot and midfoot. Thanks to the rigid outsole it guarantees the maximum stability and immobilization of the foot, optimizing thus the recovery times. It is also used as a postoperative device for all surgeries on the forefoot. Thanks to its very wide fit, the footwear can be used also with important bandages. The height of the outsole of 20 mm allows to maintain the physiological posture and to minimize the discrepancy with the contralateral limb, helping thus to simulate a natural walk. The S-Tal Stelt insole is included in the package and has an 8° talism angle, which allows to transfer the pressure from the forefoot to the hindfoot.

TECNICA CT STELT

Certification no. 2702904

Color Black

Sizes from 35 to 48

Ref. no. EX021190



Sizes						
35-36	37-38	39-40	41-42	43-44	45-46	47-48



PLANTARE S-TAL STELT

Certification no. 2702903

Sizes from 35 to 48

Ref. no. EX021250



PLANTARE S-CLOUD STELT

Certification no. 2702905

Sizes from 35 to 48

Ref. no. EX021420



DIABETIC FOOT BRACES AND POST-OPERATIVE DEVICES LINE

VARIABLE VOLUME - AMBIDEXTROUS

NO INNER
STITCHINGMEDICAL
DEVICESUITABLE FOR THE
PREVENTION OF
DIABETIC FOOTPADDED
LININGROCKING
OUTSOLE

TECNICA BIOMED

Certification no. 2702898

Color Black/White

Sizes from 35 to 46

Ref. no. EX022450



Sizes POST-OPERATIVE					
35-36	37-38	39-40	41-42	43-44	45-46

TECNICA 20 2.0

Certification no. 2702883

Color Black

Sizes from 35 to 48

Ref. no. EX021990



TECNICA 16 2.0

Certification no. 2702882

Color Black

Sizes from 35 to 48

Ref. no. EX021490



PLANTARE S-CLOUD BIOMED

Certification no. 2702906

Sizes from 35 to 46

Ref. no. EX021420



INDICATIONS

- II° post-surgical treatment

CHARACTERISTICS

Super light and very stable footwear thanks to the thermoformable stiffener, with open front, particularly suitable for deformities and injuries in the forefoot area and ideal also for physiotherapy for functional recovery. The variable internal volume of the upper, which can be completely opened, facilitates the fit of bandaged or painful feet in the phase following the brace. The stiffener is customizable by the orthopedic technician so that the fit can be adapted to the patient's foot. Non-slip outsole for maximum safety during walking. Breathable upper in washable microfiber.

Sizes POST-OPERATIVE						
35-36	37-38	39-40	41-42	43-44	45-46	47-48

INDICATIONS

- Post-surgical treatment of forefoot fractures
- Hallux valgus, hammer toes
- Bursitis
- Forefoot and interdigital ulcerations

CHARACTERISTICS

Designed in order to provide protection and relief to the forefoot. The outsole in talism release the operated area from the pressure which is shifted to the metatarsal and to the heel. Allows a precocious walking, anticipating functional recovery. Wide shaped outsole and rigid stiffener for maximum stability. Variable volume. Upper in breathable microperforated fabric washable at 40°.

Sizes POST-OPERATIVE						
35-36	37-38	39-40	41-42	43-44	45-46	47-48

DIABETIC FOOT BRACES AND POST-OPERATIVE DEVICES LINE

VARIABLE VOLUME - AMBIDEXTROUS

NO INNER
STITCHINGMEDICAL
DEVICESUITABLE FOR THE
PREVENTION OF
DIABETIC FOOTPADDED
LINING

TECNICA 8

Certification no. 2702889

Color	Grey
Sizes	from 35 to 48
Ref. no.	EX021550



INDICATIONS

- Post-surgical treatment of forefoot fractures
- Hallux valgus, hammer toes
- Bursitis
- Forefoot and interdigital ulcerations

CHARACTERISTICS

Designed in order to provide protection and relief to the forefoot. The outsole in talism release the operated area from the pressure which is shifted to the metatarsal and to the heel. Allows a precocious walking, anticipating functional recovery. Wide shaped outsole and rigid stiffener for max mum stability. Variable volume. Upper in breathable microperforated fabric w shable at 40°. With front rim in order to apply the Tecnica C8 toe-cover for thermal and hygienic protection.

TECNICA C8

Certification no. 2702551

Ref. no. EX020490



Sizes POST-OPERATIVE						
35-36*	37-38*	39-40*	41-42**	43-44**	45-46**	47-48**

Small Sizes *
Big Sizes **

TECNICA MOVIDA

Certification no. 2702892

Color	Blue
Sizes	from 35 to 46
Ref. no.	EX020970



INDICATIONS



- Post-surgical treatment

CHARACTERISTICS

Cast shoe suitable for cast bandages. Upper in breathable jeans fabric with laces and robust non-slip outsole, with a 20mm high peripheral edge for protection from water.

Sizes TECNICA MOVIDA			
35-36-37	38-39-40	41-42-43	44-45-46

INTERNAL MEASURES

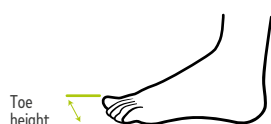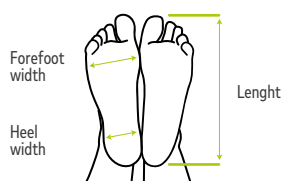


MEASURES OF THE MODELS: BIOMED, STELT, 8, 16, 20 RIGID.

Sizes	UM	35/36	37/38	39/40	41/42	43/44	45/46	47/48
Lenght	mm	235	250	260	270	285	295	310
Forefoot width	mm	85	90	95	100	105	110	115
Heel width	mm	72	74	74	78	80	82	84



PATENTED INNOVATION

SISTEMA COREX

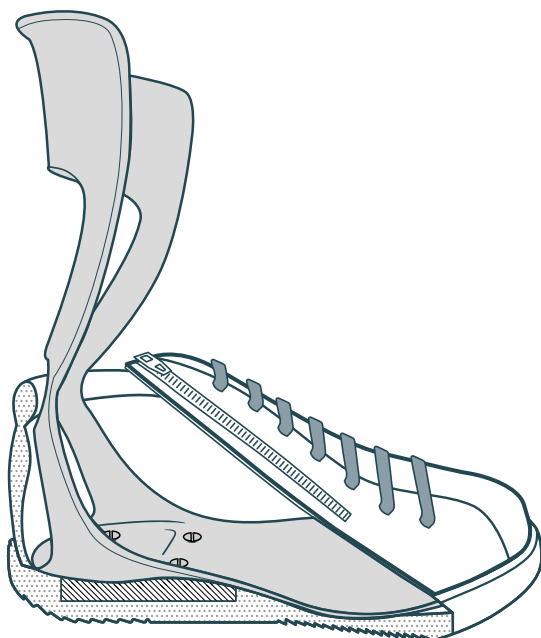

FOOTWEAR INTEGRATED WITH ORTHOPEDIC BRACE

Designed to meet the needs of patients affected by conditions such as foot drop, equinus foot, and other orthopedic conditions.

The brace and footwear become a single, integrated system, offering significant advantages in terms of comfort, stability, and performance.

WITH THE COREX SYSTEM, WE GUARANTEE...

INNOVATION IN THE SERVICE OF HEALTH: The COREX SYSTEM is a patented technology designed to ensure optimal results.

TOTAL ADAPTABILITY: The COREX SYSTEM can be applied to all products of the Energy line.

UNIVERSAL COMPATIBILITY: Compatible with all major ankle-foot orthoses (AFOs), Codivilla braces, and other types of orthopedic devices.

SUPPORT FOR PROFESSIONALS: A product designed to facilitate the work of orthopedic technicians and medical specialists.

COMPLETE KIT FOR THE ORTHOPEDIC TECHNICIAN: Includes everything needed to perform the anchoring process with precision and quality, ensuring optimal results for the patient.

WHAT MAKES THE COREX SYSTEM UNIQUE?


- 1 PATENTED ANCHORING SYSTEM
- 2 ADAPTABILITY ACROSS THE ENTIRE ENERGY LINE
- 3 USE OF CUSTOM ORTHOSES
- 4 ANCHORING KIT INCLUDED
- 5 IMPROVED PERFORMANCE
- 6 VERSATILITY

COMPLETE KIT FOR THE ORTHOPEDIC TECHNICIAN

Includes everything needed for securing the brace into the footwear

INSOLE
compensating 3 mm

INSOLE
6 mm

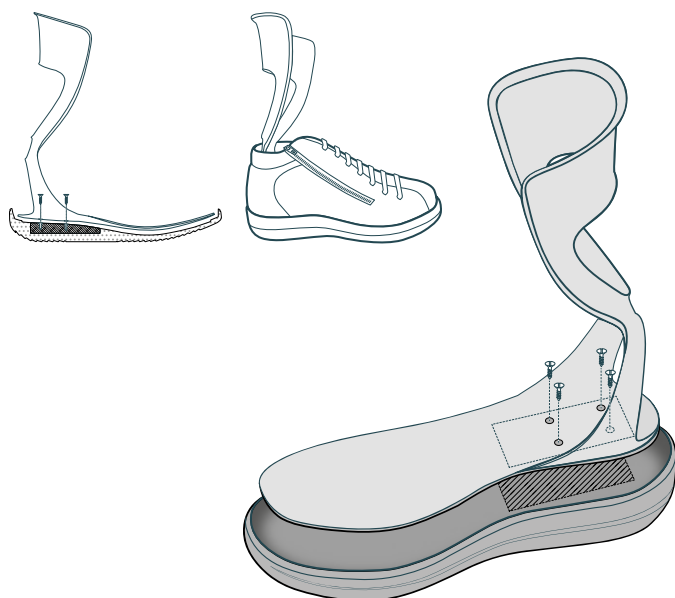

**THE FIRST
SAFETY & ORTHOPAEDIC
FOOTWEAR**



METAL FREE



SUPER LIGHT

**THE PROTECTION
FOR YOUR SENSITIVE
FEET**



EXTREME
MEDICAL WORK SOLUTION



EXTREME
MEDICAL WORK SOLUTION



INTRODUCTION

EXTREME LINE

Footwear with double certification

From the technology and experience of ExtroStyle, a leading manufacturer of footwear for foot health, comes a unique line of footwear: **EXTREME!** Especially designed to protect sensitive feet during work. **THIS FOOTWEAR POSSESSES A DOUBLE CERTIFICATION: SAFETY CERTIFICATION AS IT IS A PPE PERSONAL PROTECTIVE DEVICE ACCORDING TO EN ISO**

20345:2011 AND, AT THE SAME TIME, AN ORTHOPAEDIC FOOTWEAR MD MEDICAL DEVICE ACCORDING TO CE MDR 745/2017. They have been created for women and men who are active every day in the world of work wearing a protective device, with the opportunity to wear a shoe designed for the insertion of a customised foot orthosis.

PPE - PERSONAL PROTECTIVE EQUIPMENT



ANTI-PUNCTURE
INSERT



HEEL AREA ENERGY
ABSORPTION



ANTISTATIC



DYNAMIC
WATERPROOFING
OF THE UPPER



RESISTANCE
TO OILS AND
HYDROCARBONS



SLIPPING
RESISTANCE



200-JOULE
RESISTANT TOE CAP



BREATHABLE
UPPER

EXTREME GREY (flex outsole)

Sizes available from 36 to 48

Certification no. 2207738

REF. N° EX067250



EXTREME WHITE (flex outsole)

Sizes available from 36 to 48

Certification no. 2207903

REF. N° EX067250



The EXTREME outsole has a special tread design that guarantees high anti-slip properties, both in the longitudinal and transverse direction, as well as being antistatic and resistant to hydrocarbons, oils and greases in general.

DESCRIPTION

EXTREME LINE AS PERSONAL PROTECTIVE EQUIPMENT

EXTREME safety footwear for professional use is certified in accordance with the latest regulations; it is suitable for most types of use as it is **S3 SRC PROTECTION** grade. It is a **soft and a lightweight shoe, with leather upper and technical fabric inserts, anti-crush toe cap, antistatic and multi-layer ballistic fabric insert with the absolute guarantee of a puncture barrier, which is also highly flexible.**

TEPOR's waterproof and breathable membrane does not allow water to enter the shoe and, at the

same time, allows the foot's moisture to escape.

TE-POR thereby guarantees maximum comfort while keeping the foot dry with an optimal microclimate inside the shoe.

The **polyether sole** guarantees absolute comfort due to its lightness, ensures flexibility and absorption of vibrations and shocks, is antistatic and resistant to hydrocarbons and oils and greases in general. The extreme compactness and the special tread design guarantee excellent anti-slip properties in both longitudinal and transverse directions.

NORMATIVA EN ISO 20345:2011 CATEGORIA DI SICUREZZA		
REQUISITI MINIMI	SIGLA	S3
PUNTALE RESISTENTE A 200 JOULES	S	
RESISTENZA AGLI OLI E IDROCARBURI	FO	
ASSORBIMENTO ENERGIA NELLA ZONA TALLONE	E	
ANTISTATICA	A	
ANTIPERFORAZIONE	P	
IMPERMEABILITÀ DINAMICA DELLA TOMAIA	WRU	
RESISTENZA ALLO SCIVOLAMENTO	SR	

INDICATIONS FOR USE AS PPE

The S3 standard corresponds to work shoes with full protection, which is why they are in high demand. They combine the properties of S2 with the anti-puncture layer typical of S1 footwear. An ideal combination for many high-risk work environments, such as:

LIGHT INDUSTRY → • Food and Beverage • Pharmaceutical • Mechanical • Chemical • Automotive • Textile • Footwear • Clothing • Wood • Paper • Construction.

HEAVY INDUSTRY → • Oil and Gas • Metallurgy • Mechanical carpentry • Shipyards.

TRADES → • Personal Care • Electrician • Carpenter • Mechanic • Plumber • Mason.

LOGISTICS AND TRANSPORT → • Warehousing • Agriculture and livestock • Security and Surveillance.

CATERING AND HOSPITALITY → • Hospital Healthcare • Services and Tertiary Sector.



INDICATIONS FOR

REMOVABLE INSOLE

The standard removable insoles used for the **EXTREME** line are composed of a high-density, open-cell foam that provides total breathability while protecting against perspiration due to its moisture-absorbing and moisture-wicking qualities, keeping the foot ventilated, dry and odour-free, and providing a pleasant sensation during walking. The ergonomic design of the insole provides a high degree of comfort, minimising strain and tension in the plantar support, and cushioning impact pressure both while walking and also in a static position.

It is guaranteed that the performance has been determined by carrying out technical tests on the footwear fitted with a **removable insole**.

If it is necessary to replace the standard removable insole included with the footwear at the time of purchase, it must be replaced with an identical one supplied by the manufacturer in order to not alter the certified configuration. If it is necessary to insert a corrective and/or antalgic foot orthosis for preventive or curative foot care, in this case, the treatment shall only be performed by a qualified orthopaedic technician, who must comply with the manufacturer's recommendation that the maximum height of the orthosis near the toe cap does not exceed 6 (six) millimetres.

Compliance with the manufacturer's instructions is indispensable in order to not alter the technical tests carried out during certification.

DESCRIPTION

EXTREME S3 SRC



METAL FREE



SUPER LIGHT

The EXTREME **heel support** is padded and robust, it gently embrace the heel area while providing stability and protection.

The 1.6–1.8 cm thick **padded leather** upper, with almost no seams, accommodates the foot in a protective and shockabsorbing envelope.

The highly **breathable AIRTEK lining** provides the inside of the shoe with a complete air exchange with a "chimney effect" in order to maintain an optimal microclimate.



The **biomechanical outsole** provides a gradual, soft and continuous support from the heel to the entire sole of the foot, assuring a natural "rolling" of the foot, with vibration and shock absorption.

Outsole insert made of 100% real cork for lightness and thermal isolation.

The **ergonomic design** of the insole provides a high degree of comfort, minimising strain and tension in the plantar support, and cushioning impact pressure both while walking and also in a static position.

Composite toe cap with impact protection up to 200J and compression protection up to 15 kN.

Multi-layer **anti-puncture insert** made of ballistic fabric with the absolute guarantee of a puncture barrier, which is also highly flexible.

THE EXTREME LINE IS **METAL FREE**

The absence of metal parts does not affect the safety and the protection of the wearer; on the contrary, **it guarantees very high safety standards during any kind of performance**. EXTREME Metal Free shoes are appreciated by all those who work in environments with the presence of metal detectors, such as: airports, banks, public offices, post offices, museums, boarding areas, customs offices and much more.

DESCRIPTION

EXTREME LINE AS A MEDICAL DEVICE

Medical Device

The **EXTREME** line is especially designed to protect sensitive feet during work activities. It is produced in two versions with a FLEXIBLE outsole and with a RIGID outsole in accordance with international protocols for the treatment and care of diabetic feet. The 1.6–1.8 mm thick full-grain leather upper combined with the waterproof and breathable **TE-POR** membrane and the highly breathable **AIR-TEK** lining provide the inside of the shoe with a complete “chimney effect” air exchange while maintaining an

optimal microclimate. The absence of internal seams ensures that there are no chafing and pressure points. It envelops the anatomy of the foot and any deformities while providing comfort and preventing annoying friction. The biomechanical **EXTREME** outsole allows a natural transition between the three phases of walking: it ensures a reduction of pressure peaks on the forefoot, limits wear and tear on the joints, enables spontaneous walking and reduces the sensation of fatigue.

MD - MEDICAL DEVICE



REMOVABLE
INSOLE



MAXIMUM
SHOCK-ABSORBING
CAPACITY
EXTRA-FLEX OUTSOLE



PREDISPOSITION
FOR INSOLE



NO INTERNAL
SEAMS



MEDICAL DEVICE



SUITABLE FOR
THE PREVENTION
OF DIABETIC
FOOT DISEASE



PADDED UPPER



EXTROSTYLE S.R.L. IN ITS ROLE AS MANUFACTURER, DECLARES THAT THE EXTREME BRAND FOOTWEAR LINES ARE CLASS I NON-INVASIVE MEDICAL DEVICES, COMPLYING WITH THE REQUIREMENTS OF EU MDR 745/2017.

DIRECTIONS FOR USE

EXTREME FLEXIBLE OUTSOLE

EXTREME FLEXIBLE OUTSOLE with removable insole is a medical device used for restraint purposes with a customised insole made by qualified health professionals for the treatment of alterations in the morphology of the foot. Suitable for all those who use a customised foot orthosis to correct, prevent or treat various temporary or chronic foot disorders. Particularly indicated for the treatment of rheumatic

foot and grade 0 and 1 prevention of diabetic foot disease also in presence of neuropathies accompanied by vasculopathies, as well as metatarsalgia, polyarthritis and arthrosis. The biomechanical sole facilitates the step rolling by supporting the foot without constricting or forcing it. It reduces pressure peaks on the forefoot and limits the load on the joints.

INDICATIONS FOR USE IN THE PRESENCE OF:

- ➔ Primary prevention in diabetic foot disease with risk grade 0 or 1
- ➔ Polyarthritis
- ➔ Morton's neuroma
- ➔ Arthrosis
- ➔ Metatarsalgia





EXTREME
MEDICAL WORK SOLUTION

EXTREME LINE

INTERNAL MEASURES

FIT 12 UM mm						
SIZE	INSTEP CIRCUMFERENCE	LENGTH	FOREFOOT CIRCUMFERENCE	FOREFOOT WIDTH	HEEL WIDTH	TOE HEIGHT
48	294,5	325	296	101	71	36,5
47	290	318	291	99,5	70	36
46	285,5	311	286	98	69	35,5
45	281	304	281	96,5	68	35
44	276,5	298	277	95,5	67	34,5
43	272	291	272	94	66	34
42	267,5	284	268	93	65	32,5
41	258	278	258	89	62	32
40	293,5	271	253	87,5	61	31,5
39	249	264	249,5	86	60	31
38	244,5	258	244	84	59	30,5
37	240	251	239,5	83	57,5	30
36	235,5	244	235,5	81,5	57	29,5
35	231	238	231	80	56	29

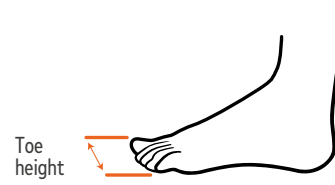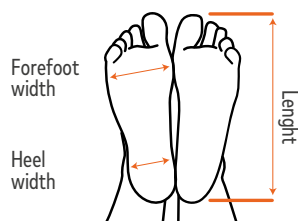

METAL FREE



SUPER LIGHT



NOTES

DATA TREATMENT

We inform you that, through the fulfillment of your order, it will be necessary to process the data you transmitted to our structure in electronic, magnetic, and paper format. The ExtroStyle company srl declares to be in compliance with the security measures provided by current legislation in the field of Privacy, Legislative Decree number 196 of the 30th June 2003 and the supplementary rules provided by the G.D.P.R (EU regulation 679/2016).

TOOLS FOR THE STORE

EXTRO STYLE PROVIDES STORES WITH A WIDE CHOICE OF DISPLAY AND SUPPORT MATERIAL



SHOE STAND



STORE DISPLAY



FOOT MEASURING DEVICE



COUNTER DISPLAYS



CATALOG



THEMATIC CATALOGS



UMBRIA Zona ind.le Colle - 06025 Nocera Umbra (PG) +39 0742 813594
MARCHE via Romagna, 37 - 62010 Morrovalle (MC) +39 0733 1825792
info@extrostyle.it www.extrostyle.it

